

ELISpot^{PLUS} for Human IgM

Product Code: 3840-2HW-Plus

CONTENTS:

Vial 1 (red top)

Monoclonal anti-IgM antibody (800 µl)

Concentration: 0.5 mg/ml

Vial 2 (yellow top)

Anti-IgM antibody, biotinylated (100 µl)

Concentration: 0.5 mg/ml

Vial 3 (white top)

Streptavidin-Horseradish Peroxidase (50 µl)

Vial 4 (black top)

Polyclonal activator: R848 (100 µl)

Concentration: 1 mg/ml

Vial 5 (blue top)

Lyophilised recombinant human IL-2 (1 µg)

To ensure total recovery of stated quantity, vials have been overfilled.

ELISpot plates (4 PVDF plates)


TMB substrate (2 x 25 ml)

STORAGE:

Shipped at ambient temperature. On arrival antibodies, streptavidin-HRP and substrate should be stored refrigerated at 4-8°C. R848 and IL-2 should be stored at -20°C. Plates may be kept at room temperature. Antibodies are supplied in sterile filtered (0.2 µm) PBS with 0.02% sodium azide. Streptavidin-HRP is supplied in PBS with 0.15% Kathon CG. R848 is supplied in sterile filtered (0.2 µm) PBS containing 2% DMSO.


Schematic illustration of the different protocols

PROTOCOL I

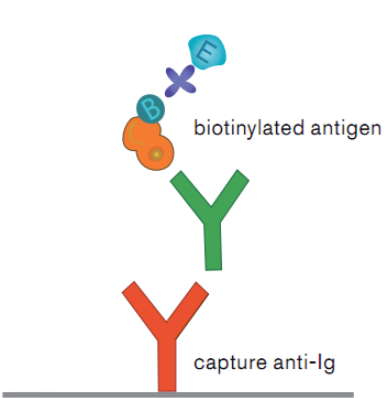

Antigen-specific Ig ELISpot

PROTOCOL II

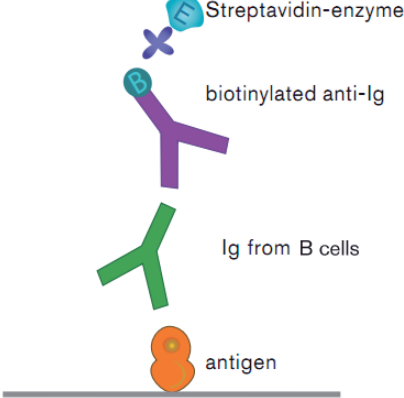

Total Ig ELISpot

Kit description

Please read through before starting the assay

This kit contains reagents for detection and enumeration of B cells secreting IgM antibodies. B cells secreting antigen-specific IgM as well as B cells secreting IgM irrespective of antigen specificity (total IgM) can be measured.

Detection of antigen-specific B cells

Detection of antigen-specific IgM-secreting B cells, visualised as spots, can be made using either of two protocols (Protocol I or II). The choice of protocol will depend on several factors such as the available amount and molecular nature of the antigen. The use of biotinylated antigen (Protocol I) may be more sensitive and result in spots of higher quality. In addition, the antigen is used in smaller quantities and is not present during cultivation.

In **Protocol I**, antibodies from all secreting B cells are captured by anti-IgM antibodies coated on the ELISpot plate. Cells producing specific antibody are then visualised by the addition of biotinylated antigen followed by Streptavidin-enzyme and substrate. Please note that antigen must be biotinylated prior to starting the experiment for this protocol.

In **Protocol II**, representing the conventional way of performing B cell ELISpot, the antigen is coated on the ELISpot plate and antigen-specific antibodies secreted by the B cells bind to the immobilised antigen. Spots are detected after addition of biotinylated anti-IgM followed by Streptavidin-enzyme and substrate.

Detection of total IgM-producing B cells

Both protocols allow determination of the total number of IgM-secreting B cells, used as a positive control or for establishing the proportion of antigen-specific B cells.

Hints and comments

Please read through before starting the assay

These suggestions are based on the detection of antigen-specific immune responses using PBMC. If using B-cell clones, mixtures of separated cell fractions etc., other protocols may have to be considered.

Plates

To obtain maximal antibody binding capacity, the plates need to first be activated by a brief treatment with ethanol. It is essential that the membrane is not allowed to dry after the treatment. If this occurs the treatment step (A2-3) needs to be repeated before adding the coating antibody.

The IgM antibodies may be polyreactive and have specificities that could cross-react with components exposed in the ELISpot wells. Hence it is important to coat the wells with irrelevant proteins or PBS alone to monitor the number of background spots.

Plate washing

Always remove the plate from the plate tray before manually emptying the plate. Washing of plates can be done using a multi-channel micropipette. In washing steps not requiring sterile conditions (C1-C5), a regular ELISA plate washer can also be used, provided that the washing head is adapted to the ELISpot plates. Avoid getting liquid on the underside of the membrane as this may cause leakage due to capillary drainage.

Cells

Both fresh and cryopreserved cells may be used with good results. However it is recommended that the latter are rested for at least one hour to allow removal of cell debris before addition to the plate. Triplicates or duplicates of 100,000-500,000 cells per well are often used to assess antigen-specific responses. For detection of total IgM, less cells (e.g. 25,000-50,000 cells per well) may be used to avoid confluent spot formation.

In vivo activated B cells

B cells that have been activated *in vivo*, for instance as a result of vaccination, may be analysed directly in the ELISpot wells without prior stimulation. Typically, cells secreting antigen-specific antibodies can be detected in the circulation 6 to 9 days after vaccination.

Memory B cells - *in vitro* activation

Memory B cells may require polyclonal stimulation before secreting detectable amounts of antibody. Pre-stimulation can be made with a mixture of R848 (1 µg/ml) and rhIL-2 (10 ng/ml) (both included in the kit) in separate tubes/plates for 24-72 hours. After pre-stimulation, wash the cells to ensure removal of any secreted antibodies. Resuspend the cells in medium and add the cell suspension to the ELISpot plate.

Reconstitution of rhIL-2: Add 1 ml PBS to obtain 1 µg/ml. Leave at room temperature for 15 minutes and then vortex. Use immediately or store in aliquots at -20°C.

Serum

The serum should be selected to support cell culture and give low background staining. We recommend the use of fetal calf serum. Alternatively serum-free medium evaluated for cell culture can be used. Human serum is not recommended as it may contain antibodies which can interfere with the assay.

Buffers

PBS for washing and dilution should be filtered (0.2 μm) for optimal results. Although possible to use, we do not recommend the inclusion of Tween or other detergents in the washing and incubation buffers.

Substrate development

Develop until distinct spots are visible in positive wells (usually 10-30 minutes). A general darkening of the membrane may occur but disappears after drying. Use water of good quality since poor quality may cause fading of TMB spots.

Biotinylation of antigen

If biotinylated antigen is to be used for detection, (Protocol I), this should be labeled before starting the ELISpot experiment. A suitable reagent for biotinylation is e.g. Sulfo-NHS-LC-Biotin (Pierce product code 21335). Further information on biotin and biotin conjugation may be found at www.piercenet.com

-The protein antigen to be biotinylated with Sulfo-NHS-LC-Biotin should be in a suitable amine-free buffer (e.g. PBS, pH 7.4). If solution contains sodium azide, this should first be removed e.g. by dialysis. For small volumes (< 200 μl) a MINI-dialysis unit Slide-A-Lyser (10.000 MWCO, Pierce) can be used. The same dialysis unit can be used for removal of free biotin (see below). Conjugation is normally performed using a 20- to 50-fold molar excess of biotin to protein where the higher ratio is for protein at low concentration.

-The Sulfo-NHS-LC-Biotin should be dissolved in pure H_2O and immediately added to the antigen solution. After incubation for 30 minutes at room temperature, free biotin may be removed by dialysis against PBS or other buffer.

The concentration of the biotinylated antigen in the assay is typically 0.01-1 $\mu\text{g}/\text{ml}$ but depends on the size of the antigen and the degree of biotinylation and should be established separately.

Protocol I

(using biotinylated antigen)

A Preparation of ELISpot plate (sterile conditions)

1. Dilute the coating anti-IgM antibody to 10 $\mu\text{g}/\text{ml}$ in sterile PBS, pH 7.4.
2. Remove the ELISpot plate (type MAIPSWU included in this kit) from the package and pre-wet the membrane by adding 50 μl 70% ethanol per well for maximum 2 minutes.
3. Wash plate 5 times with sterile water, 200 $\mu\text{l}/\text{well}$. To avoid leakage, always remove the plate from the plate tray before manually emptying the plate.
4. Add 100 $\mu\text{l}/\text{well}$ of the antibody solution and incubate overnight at 4-8°C.

B Incubation of cells in plate (sterile conditions)

See Hints and Comments section for details on cell stimulation.

1. Wash plate 5 times with sterile PBS, 200 $\mu\text{l}/\text{well}$, to remove excess antibody.
2. For blocking, add 200 $\mu\text{l}/\text{well}$ of medium containing 10% of the same serum as used for the cell suspensions. Incubate for at least 30 minutes at room temperature.
3. For analysis of *in vivo* activated B cells, remove the medium and add the cell suspension to the ELISpot plate. For analysis of memory B cells, the cells can be pre-activated separately in tubes prior to addition to the plates. See further Hints and Comments/*In vivo* activated B cells and Memory B cells.
4. Put the plate in a 37°C humidified incubator with 5% CO₂ and incubate for 16-24 hours. Do not move the plate during this time and take measures to avoid evaporation (e.g. by wrapping the plate in aluminium foil).

C Detection of spots

1. To remove the cells, empty the plate and wash 5 times with PBS, 200 $\mu\text{l}/\text{well}$.
- 2a. **Antigen-specific IgM**
Dilute biotinylated antigen to a suitable concentration (e.g. 0.01-1 $\mu\text{g}/\text{ml}$) in PBS containing 0.5% fetal calf serum (FCS). Add 100 $\mu\text{l}/\text{well}$ and incubate for 2 hours at room temperature.
- 2b. **Total IgM**
Dilute the detection antibody (anti-IgM-biotin) to 1 $\mu\text{g}/\text{ml}$ in PBS-0.5% FCS. Add 100 $\mu\text{l}/\text{well}$ and incubate for 2 hours at room temperature.
3. Wash plate as above (step C1).
4. Dilute the Streptavidin-HRP (1:1000) in PBS-0.5% FCS and add 100 $\mu\text{l}/\text{well}$. Incubate for 1 hour at room temperature.
Please note that sodium azide used in buffers will inhibit HRP activity.
5. Wash plate as above (step C1).
6. Add 100 $\mu\text{l}/\text{well}$ of ready-to-use TMB substrate solution and develop until distinct spots emerge.
7. Stop colour development by washing extensively in tap water. Remove the plate from the plate tray and rinse the underside of the membrane.
8. Leave the plate to dry. Inspect and count spots in an ELISpot reader or in a dissection microscope.
9. Store plate in the dark at room temperature.

Protocol II

(using coated antigen)

A Preparation of ELISpot plate (sterile conditions)

- | | |
|---|---|
| <p>1a. Antigen-specific IgM Dilute the antigen to suitable concentration (e.g. 1-50 µg/ml in sterile PBS, pH 7.4).</p> | <p>1b. Total IgM Dilute the coating anti-IgM antibody to 10 µg/ml in sterile PBS, pH 7.4</p> |
|---|---|
- Remove the ELISpot plate (type MAIPSWU included in this kit) from the package and pre-wet the membrane by adding 50 µl 70% ethanol per well for maximum 2 minutes.
 - Wash plate 5 times with sterile water, 200 µl/well. To avoid leakage, always remove the plate from the plate tray before manually emptying the plate.
 - Add 100 µl/well of the antigen or antibody solution and incubate overnight at 4-8°C.

B Incubation of cells in plate (sterile conditions)

See Hints and Comments section for details on cell stimulation.

- Wash plate 5 times with sterile PBS, 200 µl/well, to remove excess antibody
- For blocking, add 200 µl/well of medium containing 10% of the same serum as used for the cell suspensions. Incubate for at least 30 minutes at room temperature.
- For analysis of *in vivo* activated B cells, remove the medium and add the cell suspension to the ELISpot plate. For analysis of memory B cells, the cells can be pre-activated separately in tubes prior to addition to the plates. See further Hints and Comments/In vivo activated B cells and Memory B cells.
- Put the plate in a 37°C humidified incubator with 5% CO₂ and incubate for 16-24 hours. Do not move the plate during this time and take measures to avoid evaporation (e.g. by wrapping the plate in aluminium foil).

C Detection of spots

- To remove the cells, empty the plate and wash 5 times with PBS, 200 µl/well.
- Dilute the detection antibody (anti-IgM-biotin) to 1 µg/ml in PBS containing 0.5% fetal calf serum (PBS-0.5% FCS). Add 100 µl/well and incubate for 2 hours at room temperature.
- Wash plate as above (step C1).
- Dilute the Streptavidin-HRP (1:1000) in PBS-0.5% FCS and add 100 µl/well. Incubate for 1 hour at room temperature.
Please note that sodium azide used in buffers will inhibit HRP activity.
- Wash plate as above (step C1).
- Add 100 µl/well of ready-to-use TMB substrate solution and develop until distinct spots emerge.
- Stop colour development by washing extensively in tap water. Remove the plate from the plate tray and rinse the underside of the membrane.
- Leave the plate to dry. Inspect and count spots in an ELISpot reader or in a dissection microscope.
- Store plate in the dark at room temperature.

NOTE; for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages therefrom.



MABTECH AB
Box 1233
SE-131 28 Nacka Strand
Sweden
Tel: +46 8 716 27 00
Fax: +46 8 716 27 01
E-mail: mabtech@mabtech.com
www.mabtech.com

MABTECH Inc
M.E.B. 220
3814 West Street
Cincinnati, OH 45227
USA
Tel: +1 513 871 4500
Fax: +1 513 871 7353
E-mail: mabtech.usa@mabtech.com

MABTECH AB Büro Deutschland
Germany
Tel: +49 40 4135 7935
Fax: +49 40 4135 7945
E-mail: mabtech.de@mabtech.com

2012-02-29

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the following standards:



MABTECH AUSTRALIA Pty Ltd
44 Gresswell Road
Macleod, VIC 3085
Australia
Tel: +61 3 9459 9630
Fax: +61 3 9455 0084
E-mail: mabtech.au@mabtech.com

MABTECH AB Bureau de liaison France
BP 255, 1300 route des Crêtes
06905 Sophia Antipolis
France
Tel: +33 (0)4 92 38 80 70
Fax: +33 (0)4 92 38 80 71
E-mail: mabtech.fr@mabtech.com