

Technical Data Sheet

Canine C-Reactive Protein [CRP] ELISA Kit

Product Information

Material Number: 557826

Description

Materials Provided

Capture antibody coated microplate

One plate of 96 breakable wells coated with a goat anti-canine CRP antibody

Detection antibody/Enzyme Conjugate (100x)

120 µL of concentrated horseradish peroxidase (HRP) conjugated to a goat anti-canine CRP antibody containing stabilizers and preservative. Protect from light.

Standard (10x)

250 µL of canine serum with elevated levels of CRP.

Wash Buffer

One packet of powdered phosphate-buffered saline (PBS) with 0.05% Tween-20. Reconstitute with 1L distilled water.


TMB Substrate

A 12 mL solution containing 3,3',5,5'-tetramethylbenzidine (TMB) supplied in a protective opaque bottle. Protect from light.

Stop Solution

12 mL of diluted phosphoric acid.

The Canine C-Reactive Protein (CRP) ELISA Kit is designed for the detection and quantitation of canine CRP in canine serum. CRP is an acute-phase protein produced by the liver during conditions of inflammation, bacterial infection, or tissue trauma. Quantitation of CRP can be useful for the determination of inflammatory conditions that would be otherwise difficult to detect and monitor. This ELISA for canine CRP is a solid phase sandwich ELISA (Enzyme-Linked Immunosorbent Assay). It utilizes an antibody specific for canine CRP coated on a 96-well plate. Standards or samples are added to the wells, and any canine CRP present binds to the immobilized antibody. After the appropriate incubation, wells are washed and a horseradish peroxidase conjugated anti-canine CRP is added to produce an antibody-antigen-antibody "sandwich". Following another incubation period, the wells are washed again and a substrate solution is added, which produces a blue color in direct proportion to the amount of CRP present in the initial sample (development of a blue color indicates a positive reaction while negative reactions appear colorless). The reaction is interrupted with the Stop Solution, which changes the color from blue to yellow (negative reactions remain colorless or faintly yellow). The color is then measured at a wavelength of 450 nm (absorbance) on a spectrophotometer or on an ELISA reader.

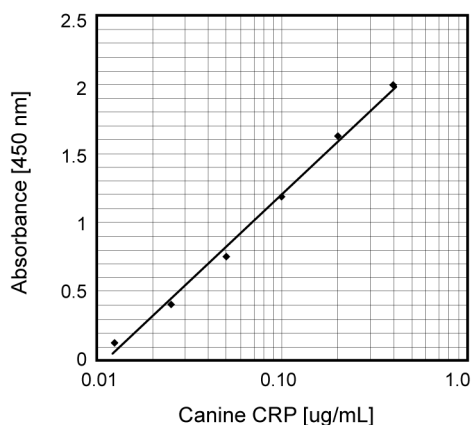

Preparation and Storage

Store undiluted at 4°C.

Kit components should be brought to room temperature (20-25°C) before opening bottles and plate pouches. Allow at least 30 minutes for this process. Do not use the kit after the expiration date.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	866.979.9408	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact


Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2015 BD

557826 Rev. 6



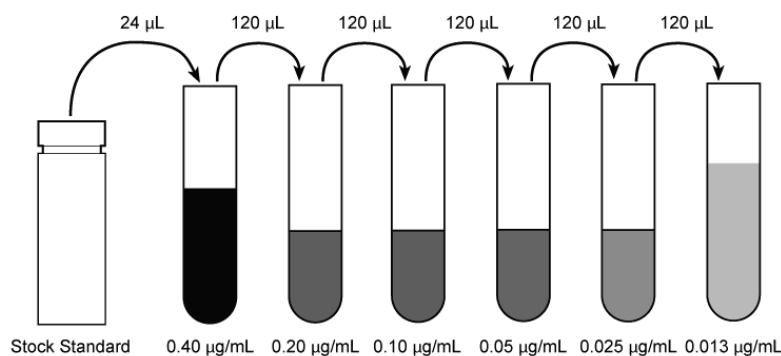

Standard #	Final concentration	Volume transferred	Diluent volume (wash buffer)	Total volume	Final volume (after dilutions)	Dilution factor
Stock (10x)	4.0 µg/mL	NA	NA	NA	NA	NA
1	0.40 µg/mL	24 µL	216 µL	240 µL	120 µL	1:10
2	0.20 µg/mL	120 µL	120 µL	240 µL	120 µL	1:2
3	0.10 µg/mL	120 µL	120 µL	240 µL	120 µL	1:2
4	0.05 µg/mL	120 µL	120 µL	240 µL	120 µL	1:2
5	0.025 µg/mL	120 µL	120 µL	240 µL	120 µL	1:2
6	0.013 µg/mL	120 µL	120 µL	240 µL	240 µL	1:2

Application Notes

Application

ELISA

Tested During Development

Recommended Assay Procedure:

Specimen Collection and Handling: Specimens should be clear, non-hemolyzed and non-lipemic. Blood samples should be collected using approved venipuncture techniques. Allow samples to clot and separate serum by centrifugation. Alternatively, use a serum separator tube (BD Vacutainer® Cat.No. 366430) and allow samples to clot for 30 minutes, then centrifuge for 10 minutes at 1000 x g. Remove serum and assay immediately or store samples at < -20°C. Avoid repeated freeze-thaw cycles.

Additional Materials Required

1. Distilled or deionized water
2. Graduated cylinder, one liter
3. Wash bottle or automated washer
4. Precision pipettes
5. Tubes or microtiter plate to prepare standard dilutions
6. Plate sealers or parafilm
7. Laboratory timer
8. Graph paper or automated data reduction software
9. Microplate reader capable of measuring absorbance at 450 nm

Recommended Assay Procedure

1. Prepare Wash Buffer by dissolving 1 packet of powdered PBS in 1L of distilled water.
2. Prepare the CRP standard by diluting the provided 10x stock standard 1:10. (Example of the dilution scheme is in the figure above). Thereafter, perform 1:2 serial dilutions five times. For example, to prepare standard for one well:
 Standard #1: Dilute the stock standard 1:10 (i.e. 24 µL of 10x stock standard mixed with 216 µL of Wash Buffer)
 Standard #2: Dilute Standard #1 1:2 (i.e. 120 µL of Standard #1 mixed with 120 µL of Wash Buffer)
 Standard #3: Dilute Standard #2 1:2 (i.e. 120 µL of Standard #2 mixed with 120 µL of Wash Buffer)
 Standard #4: Dilute Standard #3 1:2 (i.e. 120 µL of Standard #3 mixed with 120 µL of Wash Buffer)
 Standard #5: Dilute Standard #4 1:2 (i.e. 120 µL of Standard #4 mixed with 120 µL of Wash Buffer)
 Standard #6: Dilute Standard #5 1:2 (i.e. 120 µL of Standard #5 mixed with 120 µL of Wash Buffer)
3. Serum samples will likely need to be diluted before testing and a titration of the sample is strongly recommended. The suggested dilution for canine serum is 1:500 (i.e. 2 µL of serum into 1 mL Wash Buffer).
4. Add 100 µl to each well and incubate at room temperature for 30 minutes.
5. Wash the plate 3-4 times with a gentle stream of wash buffer from a wash bottle or a plate washer. Tap plates on a stack of absorbent paper towels to remove residual buffer.

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2015 BD

557826 Rev. 6



6. Prepare a 1x working concentration of the detection antibody/enzyme conjugate by diluting the provided 100x detection antibody/enzyme conjugate 1:100 with Wash Buffer (e.g. add 50 µL of 100x detection antibody/enzyme conjugate to 5 mL of Wash Buffer).
7. To each microwell being testing, add 100 µL of the 1x detection antibody/enzyme conjugate.
8. Cover the plate and incubate for 30 minutes at room temperature.
9. Wash the plate as described in step 5.
10. Add 100 µL of TMB substrate solution and incubate 5-10 minutes at room temperature.
11. Stop reaction by adding 100 µL of Stop solution.
12. Read absorbance at 450 nm within 30 minutes of stopping the reaction. If wavelength correction is available, absorbance at 630 nm may be subtracted from absorbance at 450 nm.
13. The provided 10x stock standard represents pre-diluted serum. To determine sample CRP concentrations, multiply measured values by dilution factor used (eg. by 500, if a 1:500 dilution was performed) to account for the dilution during sample preparation.

Limitations of the Procedure

1. The BD™ ELISA canine CRP kit is intended for use as an integral unit. Do not mix reagents from different kit batches.
2. Samples that generate absorbance values higher than the standard curve should be diluted and retested. Occasionally, antigen excess may be encountered in serum with high CRP values. In this situation, all the available CRP may not have reacted with the detection antibody and testing the serum at higher dilutions (e.g. 1:1000, 1:2000, or 1:4000) may be necessary.
3. Interference by drug metabolites, soluble receptors, or other binding proteins in specimens has not been thoroughly investigated. The possibility of interference cannot be excluded.
4. Do not use components past the expiration date.

Handle all serum and plasma specimens in accordance with NCCLS guidelines for preventing transmission of blood-borne infections. The control serum has not been screened for infectious agents. Since no testing can assure the complete absence of infectious agents, serum specimens and equipment coming in contact with these specimens should be handled with good laboratory practices.

Warning: PBS-T Wash Buffer Powder (component 51-915X001) contains 75.8% sodium chloride (w/w) and 1.9% (w/w) potassium chloride (w/w).

Hazard statements

May be harmful if swallowed.

Precautionary statements

Call a POSITION CENTER/doctor if you feel unwell.

Danger: Stop Solution (component 51-946P001) contains 15.23% phosphoric acid (w/w).

Hazard statements

Causes severe skin burns and eye damage.

Precautionary statements

Wear protective gloves / eye protection.

Wear protective clothing.

IF ON SKIN (or hair): Remove/Take off immediately all contaminated clothing. Rinse skin with water/shower.

IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

IF INHALED: Remove victim to fresh air and keep at rest in a position comfortable for breathing.

Dispose of contents/container in accordance with local/regional/national/international regulations.

Product Notices

1. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

References

Burton SA, Honor DJ, Mackenzie AL, Eckersall PD, Markham RJ, Horney BS. C-reactive protein concentration in dogs with inflammatory leukograms. *Am J Vet Res.* 1994; 55(5):613-618. (Biology)

Otake K, Sugimoto T, Jinbo T, et al. Physiological levels of C-reactive protein in normal canine sera. *Vet Res Commun.* 1998; 22(2):77-85. (Biology)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	866.979.9408	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2015 BD

557826 Rev. 6

