

## Dynabeads® CD25

Catalog no. 11157D

Store at 2°C to 8°C

Rev. Date: March 2012 (Rev. 003)

### Product Contents

Product contents	Volume
Dynabeads® CD25	5 mL

#### Product capacity

Whole blood: 200 mL  
MNC:  $\sim 5 \times 10^9$  cells

Dynabeads® CD25 contains  $4 \times 10^8$  beads/mL in phosphate buffered saline (PBS), pH 7.4, with 0.1% bovine serum albumin (BSA) and 0.02% sodium azide as a preservative.

**Caution:** Sodium azide may react with lead and copper plumbing to form highly explosive metal azides.

### Product Description

Isolate or deplete human CD25<sup>+</sup> T cells directly from whole blood, buffy coat or mononuclear cells (MNC) with Dynabeads® CD25. For rapid and consistent results in protein or gene expression analysis, lyse the CD25<sup>+</sup> T cells while they are still attached to the beads and directly process for further molecular analysis. The beads are mixed with the cell sample in a tube. The beads bind to the target cells during a short incubation, and then the bead-bound cells are separated by a magnet (fig. 1).

**Depletion** – Discard the bead-bound cells and use the remaining, untouched cells for any application.

**Positive isolation** – Discard the supernatant and use the bead-bound cells for downstream molecular applications.

### Downstream Applications

CD25<sup>+</sup> cells can be efficiently depleted from a sample. For rapid and consistent results in protein or gene expression analysis, lyse the CD25<sup>+</sup> T cells while still attached to the beads and directly process for further molecular analysis. For positive isolation for functional studies, cell activation/expansion, or for flow cytometer analysis, the cells need to be released after isolation. For this, we recommend using Dynabeads® FlowComp™ Flexi with your own CD25 antibody (bead-free cells). For isolation of CD4<sup>+</sup>CD25<sup>+</sup> regulatory T cells, see "Related Products".

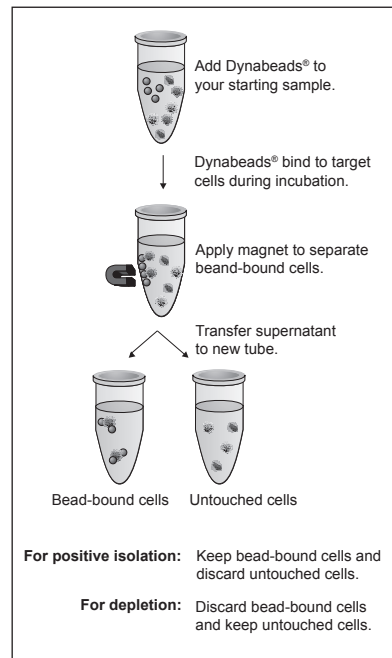


Figure 1: Overview of method

### Required Materials

- Magnet (DynaMag™ portfolio). See [www.lifetechnologies.com/magnets](http://www.lifetechnologies.com/magnets) for recommendations.
- Mixer allowing tilting and rotation of tubes (e.g. HulaMixer® Sample Mixer).
- Isolation Buffer: Ca<sup>2+</sup> and Mg<sup>2+</sup> free PBS supplemented with 0.1% BSA and 2 mM EDTA, pH 7.4.  
**Note:** BSA can be replaced by human serum albumin (HSA) or fetal calf serum (FCS). EDTA can be replaced by sodium citrate.

### General Guidelines

- Visit [www.lifetechnologies.com/samplepreparation](http://www.lifetechnologies.com/samplepreparation) for recommended sample preparation procedures.
- Use a mixer that provides tilting and rotation of the tubes to ensure that beads do not settle in the tube.
- This product should not be used with the MPC™-1 magnet (Cat. no. 12001D).
- Avoid air bubbles (foaming) during pipetting.
- Carefully follow the recommended pipetting volumes and incubation times.
- Keep all buffers cold.

### Protocol

Approximately  $1 \times 10^6$  T cells are present/mL human blood, and about 5–10% of these T cells strongly express the CD25 antigen. This protocol describes depletion or positive isolation of CD25<sup>+</sup> T cells from MNC using Dynabeads® CD25.

### Wash the Beads

See Table 1 for volume recommendations.

1. Resuspend the beads in the vial (i.e. vortex for >30 sec, or tilt and rotate for 5 min).
2. Transfer the desired volume of beads to a tube.
3. Add the same volume of Isolation Buffer, or at least 1 mL, and resuspend.
4. Place the tube in a magnet for 1 min and discard the supernatant.
5. Remove the tube from the magnet and resuspend the washed beads in the same volume of Isolation Buffer as the initial volume of beads (step 2).

### Prepare Cells

- Cells can be directly isolated from any sample such as whole blood, bone marrow, MNC suspensions or tissue digests.
- Prepare MNC to  $2.5 \times 10^7$  cells/mL in Isolation Buffer.
- See "General Guidelines" for sample preparation procedures.

## Deplete or Positively Isolate CD25<sup>+</sup> T Cells

The protocol is based on 1 mL ( $2.5 \times 10^7$ ) MNC or 1 mL whole blood/buffy coat as starting sample, but is scalable from  $2.5 \times 10^7$  –  $5 \times 10^8$  (1–50 mL). When working with lower volumes than 1 mL, use the same volumes as indicated for 1 mL. When working with larger volumes, scale up all volumes accordingly, as shown in Table 1.

1. Transfer 1 mL cells ( $2.5 \times 10^7$ ) to a tube and add 25  $\mu$ L pre-washed and re-suspended beads.
  2. Incubate for 20 min (positive isolation) or 30 min (depletion) at 2°C to 8°C with gentle tilting and rotation.
  3. Place the tube in a magnet for 2 min.
  4. For *depletion*; transfer supernatant to a new tube for further use and discard the beads.  
*or*  
For *positive isolation*; while the tube is still in the magnet, carefully remove and discard the supernatant.
  5. Remove the tube from the magnet and add 1 mL Isolation Buffer, pipet 2–3 times (or vortex 2–3 sec) and place the tube in a magnet for 2 min. While the tube is still in the magnet, carefully remove and discard the supernatant.
  6. Repeat step 5 at least once to wash the bead-bound CD25<sup>+</sup> T cells. This step is critical to obtain a high purity of isolated cells.
  7. Resuspend the cell pellet in preferred cell medium.
- Keep the cells on 2°C to 8°C until further use in downstream applications.

Table 1: Volumes for isolation/depletion of human CD25<sup>+</sup> T cells. This protocol is scalable from  $2.5 \times 10^7$  to  $5 \times 10^8$  cells.

Step	Step description	Small scale (1X)	Large scale (10X)
	Recommended tube size	5 mL	15 mL
	Recommended magnet	DynaMag™-5	DynaMag™-15
1*	Sample volume (MNC/blood/buffy)	1 mL	10 mL
1**	Bead volume	25 $\mu$ L	250 $\mu$ L
5–6	<b>For positive isolation only:</b> Wash cells (Isolation Buffer)	3 $\times$ ~1 mL	3 $\times$ ~10 mL

\*  $2.5 \times 10^7$  MNC/mL.

\*\* If very high cell-depletion efficiency is required, increase the beads volume up to double the recommended amount.

## Description of Materials

Dynabeads® CD25 are uniform, superparamagnetic polystyrene beads (4.5  $\mu$ m diameter) coated with a primary monoclonal mouse IgM antibody specific for the CD25 membrane antigen, which is predominantly expressed on human activated T cells, activated B cells and regulatory T cells.

## Related Products

Product	Cat. no.
DynaMag™-5	12303D
DynaMag™-15	12301D
DynaMag™-50	12302D
HulaMixer® Sample Mixer	15920D
Dynabeads® Regulatory CD4 <sup>+</sup> CD25 <sup>+</sup> T Cell Kit	11363D
Dynabeads® FlowComp™ Flexi	11061D

**REF** on labels is the symbol for catalog number.

## Limited Use Label License

The purchase of this product conveys to the purchaser the limited, nontransferable right to use the purchased amount of the product only to perform internal research for the sole benefit of the purchaser. No right to resell this product or any of its components is conveyed expressly, by implication, or by estoppel. This product is for internal research purposes only and is not for use in commercial applications of any kind, including, without limitation, quality control and commercial services such as reporting the results of purchaser's activities for a fee or other form of consideration. For information on obtaining additional rights, please contact [outlicensing@lifetech.com](mailto:outlicensing@lifetech.com) or Out Licensing, Life Technologies, 5791 Van Allen Way, Carlsbad, California 92008.

Manufactured by Life Technologies AS, Norway. Life Technologies AS complies with the Quality System Standards ISO 9001:2008 and ISO 13485:2003.

## Limited Product Warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at [www.lifetechnologies.com/termsandconditions](http://www.lifetechnologies.com/termsandconditions). If you have any questions, please contact Life Technologies at [www.lifetechnologies.com/support](http://www.lifetechnologies.com/support).

SPEC-06006

©2012 Life Technologies Corporation. All rights reserved. The trademarks mentioned herein are the property of Life Technologies Corporation or their respective owners, except where otherwise stated. LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

For support visit [www.lifetechnologies.com/support](http://www.lifetechnologies.com/support) or email [techsupport@lifetech.com](mailto:techsupport@lifetech.com)

[www.lifetechnologies.com](http://www.lifetechnologies.com)

 life technologies™