


## Anti-Human CD103 (Integrin alpha E) Purified

Catalog Number: 14-1038

Also Known As: Integrin  $\alpha$ E, ITGAE


RUO: For Research Use Only

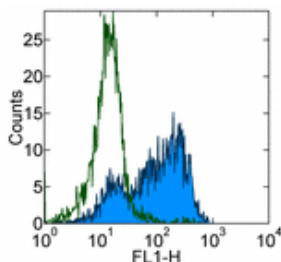

Staining of 3-day PHA-stimulated normal human peripheral blood cells with 0.25  $\mu$ g of Mouse IgG1  $\kappa$  Isotype Control Purified (cat. 14-4714) (open histogram) or 0.25  $\mu$ g of Anti-Human CD103 (Integrin  $\alpha$ E) Purified (filled histogram) followed by Anti-Mouse IgG FITC (cat. 11-4011). Total viable cells were used for analysis.

### Product Information

Contents: Anti-Human CD103 (Integrin alpha E) Purified

 Catalog Number: 14-1038

Clone: B-Ly7

Host/Isotype: Mouse IgG1,  $\kappa$

HLDA Workshop: IV B58 and V-A129

Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer



Temperature Limitation: Store at 2-8°C.



Batch Code: Refer to Vial



Use By: Refer to Vial

### Description

The B-Ly7 monoclonal antibody reacts with human CD103, the  $\alpha$ <sub>E</sub> integrin. CD103 non-covalently associates with integrin  $\beta$ <sub>7</sub>. CD103 is expressed mainly on intraepithelial lymphocytes and a small subset of peripheral lymphocytes. CD103 is also expressed by hairy cell leukemia (HCL) and by some chronic B cell lymphocytic leukemias. *In vitro* stimulation of human T cells with mitogens induces upregulation of CD103. Epithelial cell antigen, E-cadherin, binds to CD103 and mediates homing of lymphocytes to the intestinal epithelium.

### Applications Reported

B-Ly7 has been reported for use in flow cytometric analysis, and immunohistochemical staining of frozen tissue sections.

### Applications Tested

The B-Ly7 antibody has been tested by flow cytometric analysis of 2-3 day PHA-activated human peripheral blood leukocytes. This can be used at less than or equal to 0.5  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

### References

de Harven E, Soligo D, Christensen H. 1995. Immunogold labeling for the diagnosis of leukemia by transmission and scanning electron microscopy. *Scanning Microsc.* 9:1191-199

de Harven E, Christensen H, Poppema S, Scott JG. 1994. Immunogold labelling of leukemic hairy cells with the B-Ly7 monoclonal antibody: an SEM and TEM study. *Microsc Res Tech.* 28:356-67

Visser L, Dabbagh L, Poppema S. 1992. Reactivity of monoclonal antibody B-ly7 with a subset of activated T cells and T-cell lymphomas. *Hematol Pathol.* 6:37-42

Thaler J, Dietze O, Faber V, Greil R, Gastl G, Denz H, Ho AD, Huber H. 1990. Monoclonal antibody B-ly7: a sensitive marker for detection of minimal residual disease in hairy cell leukemia. *Leukemia.* 4:170-6

Schwartz R, Dienemann D, Kruschwitz M, Fritsche G, Stein H. 1990. Specificities of monoclonal antibodies B-ly7 and HML-1 are identical. *Blood.* 75:320-1.

### Related Products

14-4714 Mouse IgG1 K Isotype Control Purified

---

Not for further distribution without written consent.

Copyright © 2000-2010 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • [www.eBioscience.com](http://www.eBioscience.com) • [info@eBioscience.com](mailto:info@eBioscience.com)