


Anti-Human TNF alpha PE

Catalog Number: 12-7349

Also Known As: Tumor Necrosis Factor alpha

RUO: For Research Use Only. Not for use in diagnostic procedures.

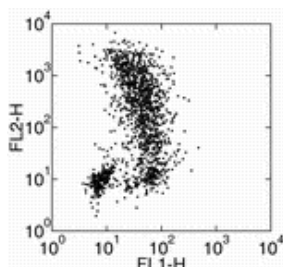

Normal human peripheral blood cells were stimulated with PMA and Ionomycin in the presence of Brefeldin A for 5 hours. The cells were surface stained with Anti-Human CD4 FITC (cat. 11-0049) and intracellularly stained with Anti-Human TNF alpha PE.

Product Information

Contents: Anti-Human TNF alpha PE


REF **Catalog Number:** 12-7349

Clone: MAb11


Concentration: ug size: 0.2 mg/ml; test size: 5 uL (0.25 ug)/test

Host/Isotype: Mouse IgG1, kappa

Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

 **Temperature Limitation:** Store at 2-8°C. Do not freeze. Light sensitive material.

LOT **Batch Code:** Refer to Vial

 **Use By:** Refer to Vial

 **Caution, contains Azide**

Description

The MAb11 antibody reacts with human tumor necrosis factor-alpha (TNF alpha), a 17 kDa cytokine produced by monocytes, macrophages, neutrophils, NK cells and CD4+ T cells. TN alpha has cytolytic activity against a range of tumor cells and is important in immune regulation. TNF alpha forms dimers and trimers and also exists as a 26 kDa membrane bound form. The MAb11 antibody is a neutralizing antibody.

Applications Reported

MAb11 is used for intracellular staining and flow cytometric analysis.

Applications Tested

This MAb11 antibody is offered in 2 formats:

- μ g size: has been tested by intracellular staining and flow cytometric analysis. This can be used at less than or equal to 0.5 μ g per test. A test is defined as the amount (μ g)/test of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

- test size: has been pre-titrated and tested by intracellular staining and flow cytometric analysis. This can be used at test size: 5 μ L (0.25 μ g)/per test The pre-titrated test size contains BSA and/or gelatin for protein stabilization.

References

Andersen H, Rossio JL, Coalter V, Poore B, Martin MP, Carrington M, Lifson JD. Characterization of rhesus macaque natural killer activity against a rhesus-derived target cell line at the single-cell level. Cell Immunol. 2004 Sep-Oct;231(1-2):85-95 (**MAb11**, crossreactivity to rhesus)

Danis VA, Franic GM, Rathjen DA, Brooks PM. Effects of granulocyte-macrophage colony-stimulating factor (GM-CSF), IL-2, interferon-gamma (IFN-gamma), tumor necrosis factor-alpha (TNF-alpha) and IL-6 on the production of immunoreactive IL-1 and TNF-alpha by human monocytes. Clin Exp Immunol. 1991 Jul;85(1):143-50.

Rathjen DA, Cowan K, Furphy LJ, Aston R. Antigenic structure of human tumour necrosis factor: recognition of distinct regions of TNF alpha by different tumour cell receptors. Mol Immunol. 1991 Jan-Feb;28(1-2):79-86.

Related Products

11-0049 Anti-Human CD4 FITC (RPA-T4)

Not for further distribution without written consent.

Copyright © 2000-2010 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • www.eBioscience.com • info@eBioscience.com