

# ELISA for Mouse IL-17A

Product Code: 3521-1H-20

---

CONTENTS, development kit for 20 plates:

**Vial 1 (yellow top)**

Monoclonal antibody IL17-I (500  $\mu$ l)

Concentration: 0.5 mg/ml

**Vial 2 (red top)**

Biotinylated monoclonal antibody IL17-II (250  $\mu$ l)

Concentration: 0.5 mg/ml

**Vial 3 (white top)**

Streptavidin-Horseradish Peroxidase (250  $\mu$ l)

**Vial 4**

Recombinant mouse IL-17A standard (1  $\mu$ g)

To ensure total recovery of stated quantity, vials have been overfilled.

**STORAGE:**

Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored at -20°C.


# General

**Intended use:** For quantitative determination of native mouse IL-17A in solution, e.g. cell culture supernatant.

**Reagents:** Antibodies are supplied in sterile-filtered (0.2  $\mu\text{m}$ ) PBS with sodium azide (0.02%). Streptavidin-HRP is supplied in PBS with 1% BSA and 0.15% Kathon CG.

**Standard range:** 5-500 pg/ml

**Standard calibration:** No international standard exists for calibration.

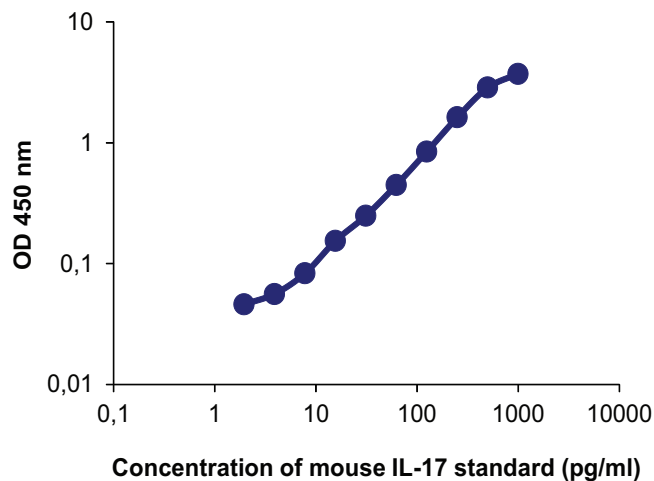

# Guidelines for Mouse IL-17A ELISA

- Day 1**
1. Coat a high protein binding ELISA plate with mAb IL17-I, diluted to 1 µg/ml in PBS, pH 7.4, by adding 100 µl/well. Incubate overnight at 4-8°C.
- Day 2**
2. Wash twice with PBS (200 µl/well).
  3. Block plate by adding 200 µl/well of PBS with 0.05% Tween 20 containing 0.1% BSA (incubation buffer). Incubate for 1 hour at room temperature.
  4. Wash five times with PBS containing 0.05% Tween.
  5. Prepare mouse IL-17A standard by reconstituting contents of vial 4 in 1 ml PBS to make a stock solution of 1 µg/ml. Allow the standard to dissolve for 5 minutes, mix thoroughly and aliquot. Store at -20°C and avoid repeated freeze-thaw cycles of the standard aliquots. For the test, prepare dilutions of the stock using the standard range as a guideline.
  6. Add 100 µl/well of samples or standards diluted in incubation buffer and incubate for 2 hours at room temperature.
  7. Wash as in step 4.
  8. Add 100 µl/well of mAb IL17-II-biotin at 0.5 µg/ml in incubation buffer or ELISA diluent for serum/plasma samples. Incubate for 1 hour at room temperature.
  9. Wash as in step 4.
  10. Add 100 µl/well of Streptavidin-HRP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature. **Please note that sodium azide used in buffers will inhibit HRP activity.**
  11. Wash as in step 4.
  12. Add 100 µl/well of appropriate substrate solution.
  13. Measure the optical density in an ELISA reader after suitable developing time. If required stop the reaction first.

**NOTE; for research use only.**

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages therefrom.



MABTECH AB  
Box 1233  
SE-131 28 Nacka Strand  
Sweden  
Tel: +46 8 716 27 00  
Fax: +46 8 716 27 01  
E-mail: [mabtech@mabtech.com](mailto:mabtech@mabtech.com)  
[www.mabtech.com](http://www.mabtech.com)

MABTECH Inc  
M.E.B. 220  
3814 West Street  
Cincinnati, OH 45227  
USA  
Tel: +1 513 871 4500  
Fax: +1 513 871 7353  
E-mail: [mabtech.usa@mabtech.com](mailto:mabtech.usa@mabtech.com)

MABTECH AB Büro Deutschland  
Germany  
Tel: +49 40 4135 7935  
Fax: +49 40 4135 7945  
E-mail: [mabtech.de@mabtech.com](mailto:mabtech.de@mabtech.com)

2013-01-24

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the following standards:



MABTECH AUSTRALIA Pty Ltd  
resolvingIMAGES  
Unit 22, 196 Settlement Road  
Thomastown Victoria 3074  
Australia  
Tel: +61 3 9466 4007  
Fax: +61 3 9466 4003  
E-mail: [mabtech.au@mabtech.com](mailto:mabtech.au@mabtech.com)

MABTECH AB Bureau de liaison France  
BP 255, 1300 route des Crêtes  
06905 Sophia Antipolis  
France  
Tel: +33 (0)4 92 38 80 70  
Fax: +33 (0)4 92 38 80 71  
E-mail: [mabtech.fr@mabtech.com](mailto:mabtech.fr@mabtech.com)