


Anti-Human CD32 PerCP-eFluor® 710

Catalog Number: 46-0329

Also known as: Low Affinity IgG Receptor 2a, FcGR2a, IGFR2

RUO: For Research Use Only. Not for use in diagnostic procedures.

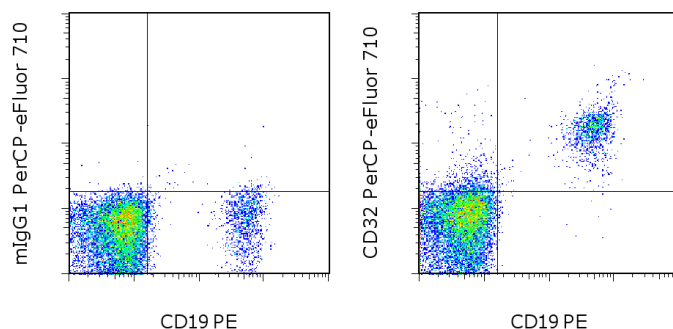

Staining of normal human peripheral blood cells with Anti-Human CD19 PE (cat. 12-0199) and Mouse IgG1 K Isotype Control PerCP-eFluor® 710 (cat. 46-4714) (left) or Anti-Human CD32 PerCP-eFluor® 710 (right). Cells in the lymphocyte gate were used for analysis.

Product Information

Contents: Anti-Human CD32 PerCP-eFluor® 710



Catalog Number: 46-0329

Clone: 6C4 (CD32)

Concentration: 5 uL (0.5 ug)/test

Host/Isotype: Mouse IgG1, kappa



Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

Temperature Limitation: Store at 2-8°C. Do not freeze. Light-sensitive material.

Batch Code: Refer to vial

Use By: Refer to vial

Contains sodium azide

Description

This 6C4 monoclonal antibody reacts with human CD32 (also known as Fc gammaRII). This antibody recognizes two isoforms of the receptor, Fc gammaRIIA and Fc gammaRIIB. These 40-kDa polymorphic transmembrane glycoproteins are expressed on B cells, granulocytes, monocytes, macrophages, and platelets. Moreover, these receptors are detected on natural killer cells. CD32 enables interaction between Fc gammaRII-expressing cells and opsonized antigen or IgG-containing immune complexes. Both receptors exhibit low affinity towards IgG and play a role in inflammation and autoimmune disease. This clone has been reported to inhibit Ig binding in a rosette blocking assay.

Applications Reported

This 6C4 (CD32) antibody has been reported for use in flow cytometric analysis.

Applications Tested

This 6C4 (CD32) antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (0.5 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.

PerCP-eFluor® 710 can be used in place of PE-Cy5, PE-Cy5.5 or PerCP-Cy5.5. PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm). Please make sure that your instrument is capable of detecting this fluorochrome. For a filter configuration, we recommend using the 685 LP dichroic mirror and 710/40 band pass filter, however the 695/40 band pass filter is an acceptable alternative.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 uL cell sample + 100 uL IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • www.ebioscience.com • info@ebioscience.com

Anti-Human CD32 PerCP-eFluor® 710

Catalog Number: 46-0329

Also known as: Low Affinity IgG Receptor 2a, FcGR2a, IGFR2

RUO: For Research Use Only. Not for use in diagnostic procedures.

References

Fossati G, Bucknall RC, and Edwards SW. Fc gamma receptors in autoimmune diseases. Eur. J. Clin. Invest. 2001 31(9):821-831.

Vely F, Gruel N, Moncuit J, Cochet O, Rouard H, Dare S, Galon J, Sautes C, Fridman WH, Teillaud JL. A new set of monoclonal antibodies against human Fc gamma RII (CD32) and Fc gamma RIII (CD16): characterization and use in various assays. Hybridoma 1997 Dec; 16(6):519-28. (6C4, FC, FA (blocking), Pubmed)

Related Products

12-0199 Anti-Human CD19 PE (HIB19)

46-4714 Mouse IgG1 K Isotype Control PerCP-eFluor® 710 (P3.6.2.8.1)