


## Anti-Mouse CD3e eFluor® 450

**Catalog Number:** 48-0033

**Also known as:** CD3 epsilon

**RUO: For Research Use Only. Not for use in diagnostic procedures.**

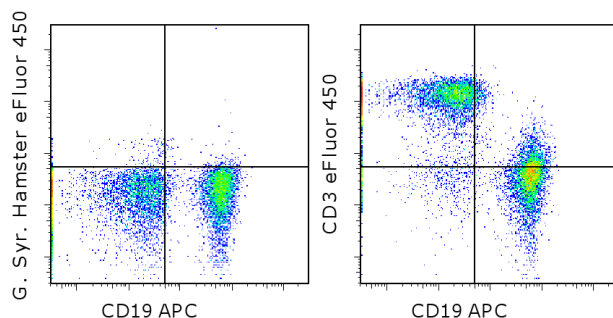

Staining of C57BL/6 splenocytes with Anti-Mouse CD19 PE (cat. 12-0193) and 0.5 ug of Anti-Mouse CD3e eFluor® 450. Total viable cells were used for analysis.

### Product Information



**Contents:** Anti-Mouse CD3e eFluor® 450

**Catalog Number:** 48-0033

**Clone:** eBio500A2 (500A2)

**Concentration:** 0.2 mg/mL

**Host/Isotype:** Golden Syrian Hamster IgG



**Formulation:** aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

**Temperature Limitation:** Store at 2-8°C. Do not freeze. Light sensitive material.



**Batch Code:** Refer to vial



**Use By:** Refer to vial

### Description

The eBio500A2 monoclonal antibody reacts with the 25 kD epsilon subunit of the mouse CD3 complex. CD3 subunits  $\gamma$ ,  $\delta$  and  $\epsilon$  are required for proper assembly, trafficking and surface expression of the TCR complex. CD3 is expressed by thymocytes in a developmentally regulated manner and by all mature T cells, as well as NKT cells. Binding of eBio500A2 to CD3 initiates the intracellular biochemical pathway resulting in cellular activation and proliferation. The 500A2 antibody is able to partially cross-block the 17A2 and 145-2C11 antibodies, indicating that all three of these anti-CD3 antibodies recognize distinct but overlapping epitopes.

### Applications Reported

This eBio500A2 (500A2) antibody has been reported for use in flow cytometric analysis.

### Applications Tested

This eBio500A2 (500A2) antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 1  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from  $10^5$  to  $10^8$  cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

**eFluor® 450 is a replacement for Pacific Blue®. eFluor® 450 emits at 456 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochrome.**

### References

Havran WL, Poenie M, Kimura J, Tsien R, Weiss A, Allison JP. Expression and function of the CD3-antigen receptor on murine CD4+8+ thymocytes. *Nature*. 1987 Nov 12-18;330(6144):170-3. (500A2, FC, FA, PubMed)

Portoles P, Rojo J, Golby A, Bonneville M, Gromkowski S, Greenbaum L, Janeway CA Jr, Murphy DB, Bottomly K. Monoclonal antibodies to murine CD3 epsilon define distinct epitopes, one of which may interact with CD4 during T cell activation. *J Immunol*. 1989 Jun 15;142(12):4169-75. (PubMed)

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • [www.ebioscience.com](http://www.ebioscience.com) •  
[info@ebioscience.com](mailto:info@ebioscience.com)

---

## Anti-Mouse CD3e eFluor® 450

**Catalog Number:** 48-0033

**Also known as:** CD3 epsilon

**RUO: For Research Use Only. Not for use in diagnostic procedures.**

---

Ortaldo JR, Winkler-Pickett R, Mason AT, Mason LH. The Ly-49 family: regulation of cytotoxicity and cytokine production in murine CD3+ cells. J Immunol. 1998 Feb 1;160(3):1158-65. (**500A2**, FC, PubMed)

### Related Products

12-0193 Anti-Mouse CD19 PE (eBio1D3 (1D3))