
Technical Data Sheet PerCP-Cy™5.5 Mouse Anti-Human CD14

Product Information	
Material Number:	550787
Size:	0.1 mg
Concentration:	0.2 mg/ml
Clone:	M5E2
Isotype:	Mouse IgG2a, ĸ
Reactivity:	QC Testing: Human Tested in Development: Baboon, Rhesus, Cynomolgus, Dog
Workshop:	III 329
Storage Buffer:	Aqueous buffered solution containing $\leq 0.09\%$ sodium azide.

Description

Reacts with the 53-55 kDa glycosylphosphatidylinositol (GPI)-anchored single chain glycoprotein expressed at high levels on monocytes. Additionally, CD14 antibody reacts with interfollicular macrophages, reticular dendritic cells, and some Langerhans cells. CD14 has been identified as a high affinity cell-surface receptor for complexes of lipopolysaccharide (LPS) and serum LPS-binding protein, LPB.

This antibody is routinely tested by flow cytometric analysis. Other applications were tested at BD Biosciences Pharmingen during antibody development only or reported in the literature.

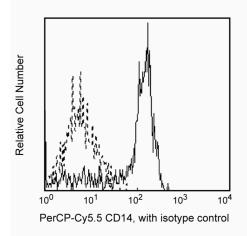

Profile of peripheral blood monocytes analyzed by flow cvtometrv

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with PerCP-Cy5.5 under optimum conditions, and unconjugated antibody and free PerCP-Cy5.5 were removed by gel filtration chromatography. Storage of PerCP-Cy5.5 conjugates in unoptimized diluent is not recommended and may result in loss of signal intensity.

Store undiluted at 4° C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application	
Flow cytometry	Routinely Tested
Suggested Companion Products	

Suggested Companion Products

Catalog Number	Name			Size	Clone	
550927	PerCP-Cy5.5 Mouse IgG2a, κ Isotype Control			0.1 mg	G155-178	
BD Biosciences						
www.bdbiosciences.com						
	32.53.720.550 0	apan)120.8555.90 vw.bdbioscience	Asia Pacific 65.6861.0633 es.com/how_to	Latin America/Caribbean 55.11.5185.9995 order/		BD
Conditions: The information disclosed of any patents. BD Biosciences will no use of our products. Purchase does no product or as a component of anothe- written authorization of Becton Dicki For Research Use Only. Not for use in BD, BD Logo and all other trademarks	t be held responsible of include or carry any r product. Any use of nson and Company is diagnostic or therape	for patent infringe y right to resell or t f this product othe s strictly prohibited eutic procedures. No	ement or other vio ransfer this produc r than the permitte ot for resale.	lations that may occur with the t either as a stand-alone ed use without the express		

Product Notices

- 1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 3. PerCP-Cy5.5–labelled antibodies can be used with FITC- and R-PE–labelled reagents in single-laser flow cytometers with no significant spectral overlap of PerCP-Cy5.5, FITC, and R-PE fluorescence.
- 4. PerCP-Cy5.5 is optimized for use with a single argon ion laser emitting 488-nm light. Because of the broad absorption spectrum of the tandem fluorochrome, extra care must be taken when using dual-laser cytometers, which may directly excite both PerCP and Cy5.5™. We recommend the use of cross-beam compensation during data acquisition or software compensation during data analysis.
- 5. CyTM is a trademark of Amersham Biosciences Limited.
- 6. This product is subject to proprietary rights of Amersham Biosciences Corp. and Carnegie Mellon University and made and sold under license from Amersham Biosciences Corp. This product is licensed for sale only for research. It is not licensed for any other use. If you require a commercial license to use this product and do not have one return this material, unopened to BD Biosciences, 10975 Torreyana Rd, San Diego, CA 92121 and any money paid for the material will be refunded.
- 7. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/pharmingen/colors.
- 8. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

References

McMichael AJ, Beverly PCL, Gilks W, et al, ed. *Leukocyte Typing III: White Cell Differentiation Antigens*. New York: Oxford University Press; 1987.(Clone-specific) Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. New York: Oxford University Press; 1995.(Biology) Knapp W, Dorken B, et al, ed. *Leucocyte Typing IV*. New York: Oxford University Press; 1989.(Biology)

Wright SD, Ramos RA, Tobias PS, Ulevitch RJ, Mathison JC. CD14, a receptor for complexes of lipopolysaccharide (LPS) and LPS binding protein. *Science*. 1990; 249(4975):1431-1433.(Biology)

BD Biosciences

 www.bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 888.259.0187
 32.53.720.550
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

 For country-specific contact information, visit
 www.bdbiosciences.com/how_to_order/
 Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product on a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

 For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
 BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2007 BD

