## Technical Data Sheet

# APC Mouse IgG2a, κ Isotype Control

#### **Product Information**

**Material Number:** 551414 Size: 50 tests 20 µl Vol. per Test: G155-178 Clone:


Immunogen: TNP-keyhole limpet hemocyanin Isotype: Mouse (BALB/c) IgG2a, κ Reactivity: Unknown specificity

QC Testing: Rhesus or Baboon or Cynomolgus

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

#### Description

The G155-178 clone has an unknown specificity. Trinitrophenol (TNP), the immunogen, is a hapten not expressed on human, mouse, rat or non-human primate cells. In the absence of specific binding, this antibody may bind non-specifically to immunoglobulin Fc receptors. The immunoglobulin secreted by the G155-178 hybridoma was selected as an mouse IgG2a, κ isotype control following screening for low background binding on a variety of mouse and human tissues.

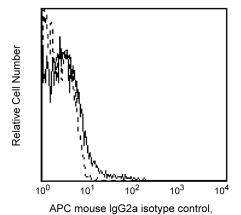

Profile of peripheral blood lymphocytes of Rhesus macaque (macaca mulatta) analyzed by flow cytometry.

# with autofluorescence

### **Preparation and Storage**

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

#### **Application Notes**

## Application

Flow cytometry	Routinely Tested
Isotype control	Routinely Tested

#### **Recommended Assay Procedure:**

An isotype control should be used at the same concentration as the antibody of interest.

### **Product Notices**

- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- Source of all serum proteins is from USDA inspected abattoirs located in the United States.

## **BD Biosciences**

bdbiosciences.com

Asia Pacific Europe 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how\_to\_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation cof any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD

