## **Technical Data Sheet**


# РЕ-Су™5 Mouse IgG2a, к Isotype Control

#### **Product Information**

Material Number:	555575
Size:	100 tests
Vol. per Test:	20 µl
Clone:	G155-178
Immunogen:	TNP-keyhole limpet hemocyanin
Isotype:	Mouse (BALB/c) IgG2a, κ
Storage Buffer:	Aqueous buffered solution containing BSA and $\leq 0.09\%$ sodium azide.

#### Description

The G155-178 clone has an unknown specificity. Trinitrophenol (TNP), the immunogen, is a hapten not expressed on human, mouse, rat or non-human primate cells. In the absence of specific binding, this antibody may bind non-specifically to immunoglobulin Fc receptors. The immunoglobulin secreted by the G155-178 hybridoma was selected as an mouse IgG2a,  $\kappa$  isotype control following screening for low background binding on a variety of mouse and human tissues.

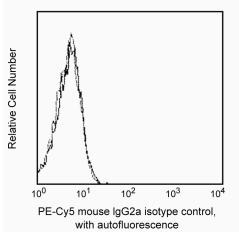

Profile of peripheral blood lymphocytes analyzed on a FACScan (BDIS, San Jose, CA)

#### **Preparation and Storage**

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with PE-Cy5 (formerly known as BD Cy-Chrome™) under optimum conditions, and unconjugated antibody and free PE-Cy5 were removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

#### Application Notes

Application							
	Flow cytometry	Routinely Tested					
	Isotype control	Routinely Tested					

#### **Product Notices**

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^{6}$  cells in a 100-µl experimental sample (a test).
- 2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 3. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 5. PE-Cy5 is optimized for use with a single argon ion laser emitting 488-nm light. Because of the broad absorption spectrum of the PE-Cy5 tandem fluorochrome, extra care must be taken when using dual-laser cytometers which may directly excite both PE and Cy5<sup>TM</sup>.
- 6. This product is subject to proprietary rights of Amersham Biosciences Corp. and Carnegie Mellon University and made and sold under license from Amersham Biosciences Corp. This product is licensed for sale only for research. It is not licensed for any other use. If you require a commercial license to use this product and do not have one return this material, unopened to BD Biosciences, 10975 Torreyana Rd, San Diego, CA 92121 and any money paid for the material will be refunded.

### **BD Biosciences**

bdbiosciences.com									
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean				
877.232.8995	888.268.5430	32.53.720.550	0120.8555.90	65.6861.0633	0800.771.7157				
For country-specific contact information, visit <b>bdbiosciences.com/how_to_order</b> /									
Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to result or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale. BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD									



- 7. Cy is a trademark of Amersham Biosciences Limited. This conjugated product is sold under license to the following patents: US Patent Nos. 5,486,616; 5,569,587; 5,569,766; 5,627,027.
- 8. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 9. Source of all serum proteins is from USDA inspected abattoirs located in the United States.