

Technical Data Sheet

FITC Mouse Anti-Human CD55


Product Information

Material Number:	555693
Alternate Name:	DAF
Size:	100 tests
Vol. per Test:	20 µl
Clone:	IA10
Isotype:	Mouse IgG2a κ
Reactivity:	QC Testing: Human
Workshop:	V BP352, S031
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.


Description

Reacts with CD55, decay-accelerating factor (DAF), a glycosylphosphatidylinositol (GPI)-anchored single chain glycoprotein of approximately 70 kDa, expressed on hematopoietic cells. It has been suggested that the role of DAF is to protect cells from damage by autologous complement preventing the amplification steps of the complement cascade by interfering with the assembly of the C3-convertases, C4b2a and C3bBb, and the C5-convertases, C4b2a3b and C3bBb3b.

This antibody is routinely tested by flow cytometric analysis. Other applications were tested at BD Biosciences Pharmingen during antibody development only or reported in the literature.

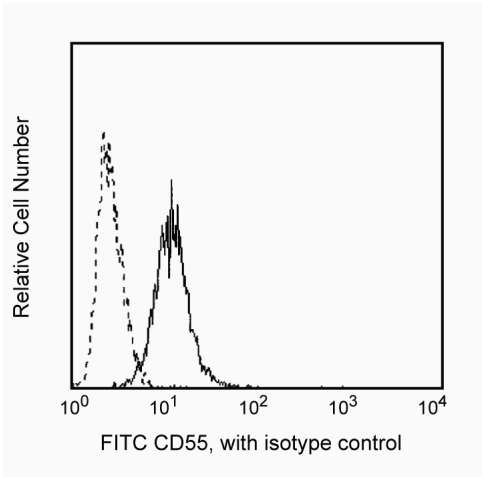

Profile of K562 cells expressing glycosylphosphatidylinositol (GPI) anchored protein analyzed on a FACScan (BDIS, San Jose, CA)

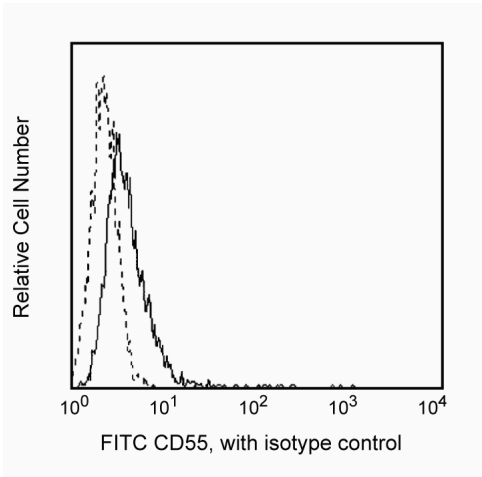

Profile of glycosylphosphatidylinositol (GPI) anchor-defective mutant cell line analyzed on a FACScan (BDIS, San Jose, CA)

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed. Store undiluted at 4° C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

Flow cytometry	Routinely Tested
----------------	------------------

BD Biosciences
bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2006 BD



Suggested Companion Products

Catalog Number	Name	Size	Clone
555573	FITC Mouse IgG2a Kappa Isotype Control	100 tests	G155-178

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 X 10⁶ cells in a 100-μl experimental sample (a test).
2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
3. Please refer to www.bdbiosciences.com/pharmlngen/protocols for technical protocols.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/pharmlngen/colors.
5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
6. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

References

Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. New York: Oxford University Press; 1995.(Clone-specific)
Kinoshita T, Medof ME, Silber R, Nussenzweig V. Distribution of decay-accelerating factor in the peripheral blood of normal individuals and patients with paroxysmal nocturnal hemoglobinuria. *J Exp Med*. 1985; 162(1):75-92.(Biology)
Loveland BE, Szokolai K, Johnstone RW, McKenzie IF. Coordinate functions of multiple complement regulating molecules, CD46, CD55, and CD59. *Transplant Proc*. 1994; 26(3):1070-1071.(Biology)