# **Technical Data Sheet**

# PE Mouse Anti-Human CD18

#### **Product Information**

**Material Number:** 555924 Alternate Name: Integrin β2 chain 100 tests Size Vol. per Test: 20 µl Clone: 6.7

**Isotype:** Mouse IgG1, κ Reactivity: QC Testing: Human


Tested in Development: Rhesus, Cynomolgus, Pig, Rabbit

Workshop:

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

# Description

Reacts with β2 integrin, a 90-95 kDa transmembrane protein, which associates non-covalently as a heterodimer with the a chains of CD11a, CD11b, or CD11c. CD18 is expressed on lymphocytes, monocytes, and more weakly on granulocytes. LFA-1 has been shown to play an important role in homotypic and heterotypic cellular adhesion in immune and inflammatory responses.

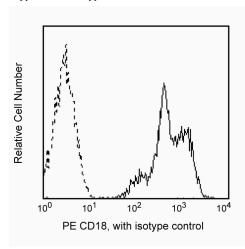

Profile of peripheral blood lymphocytes analyzed on a BD FACScan™ (BDIS, San Jose, CA)

# **Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed.

# **Application Notes**

### Application

| -FP            |                  |  |
|----------------|------------------|--|
| Flow cytometry | Routinely Tested |  |

# Suggested Companion Products

| Catalog Number | Name                             | Size      | Clone   |
|----------------|----------------------------------|-----------|---------|
| 554656         | Stain Buffer (FBS)               | 500 ml    | (none)  |
| 555749         | PE Mouse IgG1, κ Isotype Control | 100 tests | MOPC-21 |

#### **Product Notices**

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10<sup>6</sup> cells in a 100-μl experimental
- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

## **BD Biosciences**

bdbiosciences.com


Asia Pacific Europe Japan 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how\_to\_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD

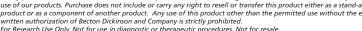



- 5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 6. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 7. An isotype control should be used at the same concentration as the antibody of interest.

#### References

Barclay NA, Brown MH, Birkeland ML, et al, ed. *The Leukocyte Antigen FactsBook*. San Diego, CA: Academic Press; 1997. (Biology) Schlossman S, Boumell L, et al, ed. *Leucocyte Typing V*. New York: Oxford University Press; 1995. (Clone-specific)

555924 Rev. 3 Page 2 of 2