

Technical Data Sheet

FITC Mouse Anti-Human Ki-67 Set**Product Information**


Material Number: 556026
Size: 100 tests
Reactivity: QC Testing: Human
 Tested in Development: Mouse, Pig, Rat

Component: 51-36524X
Description: FITC Mouse Anti-Human Ki-67
Size: 100 tests (1 ea)
Vol. per Test: 20 µl
Clone Name: B56
Immunogen: Human Ki-67
Isotype: Mouse IgG1, κ
Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.


Component: 51-35404X
Description: FITC Mouse IgG1, κ Isotype Control
Size: 100 tests (1 ea)
Vol. per Test: 20 µl
Clone Name: MOPC-21
Isotype: Mouse IgG1, κ
Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Recognizes Ki-67, a nuclear cell proliferation-associated antigen expressed in all active stages of the cell cycle. Ki-67 is revealed as a double band (345 and 395 kDa) in western blot analysis of proliferating cells. B56 was developed using an immunogen composed of the immunodominant epitope of the Ki-67 protein. Antibodies B56 and MIB 1 react with this immunogen. Flow cytometric analysis reveals that the binding of B56-PE can be blocked by MIB 1 purified antibody.

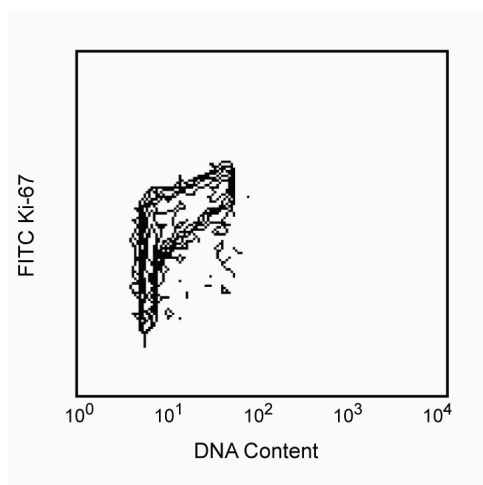

Profile of Ki-67 on MOLT-4 cells analyzed on a FACSscan (BDIS, San Jose, CA)

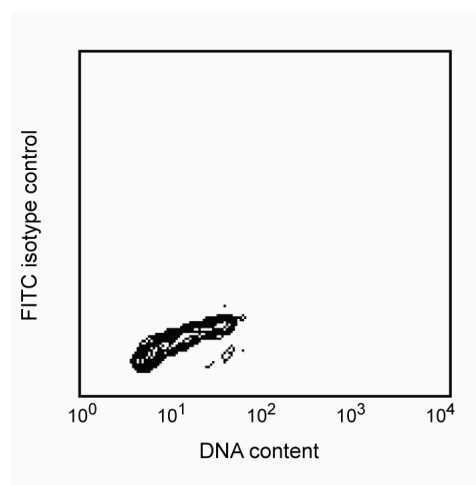

Profile of mouse IgG1 isotype control on MOLT-4 cells analyzed on a FACSscan (BDIS, San Jose, CA)

BD Biosciences

www.bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit www.bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2007 BD



Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed.

Store undiluted at 4° C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

Recommended Assay Procedure:

Ki-67 staining protocol by flow cytometry:

1. Harvest, count and pellet cells following standard procedures (Note: Ki-67 is expressed by the proliferative cells. You may get no staining with the resting cells, e.g. unstimulated PBMC).
2. While vortexing, add 5 ml drop by drop of cold 70% - 80% ethanol into the cells pellet (1-5 x 10⁶ cells). Then incubate at -20°C for 2 hours minimum. These fixed cells can be used up to 60 days after fixing (Store at -20°C).
3. Add 30-40 ml wash buffer (PBS with 1% FBS, 0.09% NaN₃ pH7.2) to the fixed cells. Centrifuge the cells for 10 minutes at 1000 rpm and aspirate supernatant. Wash one more time with 30-40 ml of wash buffer. Centrifuge at 1000 rpm for 10 minutes and aspirate supernatant.
4. Resuspend the cells to a concentration of 1 X 10⁶/ ml (1 X 10⁶/100 µl).
5. Transfer 100 µl cell suspension into each fresh tube.
6. Add 20 µl of properly diluted antibody according to the protocol into the tubes above. Mix gently.
7. Incubate the tubes at room temperature (RT) for 20-30 minutes in the dark.
8. Wash with 2 ml of PBS washing buffer at 1000 rpm for 5 minutes.
9. Aspirate the supernatant.
10. For direct conjugated antibody: go to steps 13 & 14.
11. For purified antibody: add 50 µl of diluted secondary antibody at optimal concentration (Cat. No. 555988), incubate at RT for 30 minutes in the dark.
12. Repeat step 8 & 9.
13. Add 0.5 ml of PBS wash buffer into each tube. For FITC conjugated antibody, add 10 µl of PI Staining Solution (Cat. No. 556463); for PE conjugated antibody, add 20 µl BD Via-Probe™ Cell Viability Solution (Cat. No. 555816) into each tube.
14. Analyze the sample with FACS.

Product Notices

1. This antibody has been optimized and preassayed with its matched isotype control to be used at the recommended volume of 20 µl/test. Titration of the reagents or substituting with other (non-matched) isotype control is NOT recommended.
2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

References

- Kubbutat MH, Key G, Duchrow M, Schluter C, Flad HD, Gerdes J. Epitope analysis of antibodies recognising the cell proliferation associated nuclear antigen previously defined by the antibody Ki-67 (Ki-67 protein). *J Clin Pathol.* 1994; 47(6):524-528.(Biology)
- Shi SR, Key ME, Kalra KL. Antigen retrieval in formalin-fixed, paraffin-embedded tissues: an enhancement method for immunohistochemical staining based on microwave oven heating of tissue sections. *J Histochem Cytochem.* 1991; 39(6):741-748.(Biology)