

## Technical Data Sheet


## PE Mouse Anti- BrdU Set

## Product Information

|                         |                                                                                 |
|-------------------------|---------------------------------------------------------------------------------|
| <b>Material Number:</b> | <b>556029</b>                                                                   |
| <b>Size:</b>            | 100 tests                                                                       |
| <b>Reactivity:</b>      | QC Testing: Human                                                               |
| <b>Component:</b>       | <b>51-35405X-2</b>                                                              |
| <b>Description:</b>     | PE Mouse IgG1 $\kappa$ Isotype Control                                          |
| <b>Size:</b>            | 100 tests (1 ea)                                                                |
| <b>Vol. per Test:</b>   | 20 $\mu$ l                                                                      |
| <b>Clone Name:</b>      | MOPC-21                                                                         |
| <b>Isotype:</b>         | Mouse IgG1, $\kappa$                                                            |
| <b>Storage Buffer:</b>  | Aqueous buffered solution containing BSA and $\leq 0.09\%$ sodium azide.        |
| <b>Component:</b>       | <b>51-33285X</b>                                                                |
| <b>Description:</b>     | PE Mouse Anti-BrdU (5-bromo-2'-deoxyuridine, 5-Bromouracil deoxyriboside, BUdR) |
| <b>Size:</b>            | 100 tests (1 ea)                                                                |
| <b>Vol. per Test:</b>   | 20 $\mu$ l                                                                      |
| <b>Clone Name:</b>      | 3D4                                                                             |
| <b>Isotype:</b>         | Mouse IgG1, $\kappa$                                                            |
| <b>Storage Buffer:</b>  | Aqueous buffered solution containing BSA and $\leq 0.09\%$ sodium azide.        |

## Description

Reacts with bromodeoxyuridine (BrdU) incorporated in place of thymidine during the S-phase of cycling cells. This technique enables measurement of the frequency of DNA-synthesizing cells in studies of cell activation and proliferation.

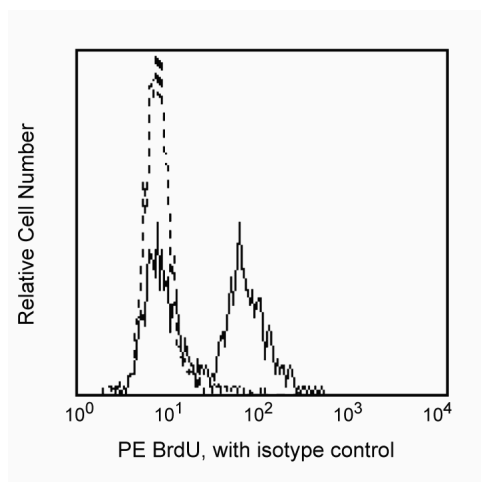

Profile of permeabilized REH cells analyzed by flow cytometry

## Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

## Application Notes

## Application

## BD Biosciences

bdbiosciences.com

| United States | Canada       | Europe        | Japan        | Asia Pacific | Latin America/Caribbean |
|---------------|--------------|---------------|--------------|--------------|-------------------------|
| 877.232.8995  | 888.259.0187 | 32.53.720.550 | 0120.8555.90 | 65.6861.0633 | 55.11.5185.9995         |

For country-specific contact information, visit [bdbiosciences.com/how\\_to\\_order/](http://bdbiosciences.com/how_to_order/)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



|                                         |                  |
|-----------------------------------------|------------------|
| Flow cytometry                          | Routinely Tested |
| Intracellular staining (flow cytometry) | Routinely Tested |

#### Recommended Assay Procedure:

Note: this product is routinely tested on fixed and permeabilized REH cell line with 70-80% cold ethanol. Please see the detail protocol at [http://www.bdbiosciences.com/pharmingen/protocols/BrdU\\_Incorporation.shtml](http://www.bdbiosciences.com/pharmingen/protocols/BrdU_Incorporation.shtml) This product is NOT recommended to use with the BrdU flow kits, e.g. Cat. No. 559619.

#### Product Notices

1. This antibody has been optimized and preassayed with its matched isotype control to be used at the recommended volume of 20 ul/test. Titration of the reagents or substituting with other (non-matched) isotype control is NOT recommended.
2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
3. Please refer to [www.bdbiosciences.com/pharmingen/protocols](http://www.bdbiosciences.com/pharmingen/protocols) for technical protocols.
4. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

#### References

Keren DF, Hanson CA, Hurtubise PE, ed. *Flow Cytometry and Clinical Diagnosis*. Chicago: American Society of Clinical Pathologists Press; 1994:1-676.(Biology)  
Dolbeare F, Gratzner H, Pallavicini MG, Gray JW. Flow cytometric measurement of total DNA content and incorporated bromodeoxyuridine. *Proc Natl Acad Sci U S A*. 1983; 80(18):5573-5577.(Biology)