

Technical Data Sheet

Alexa Fluor® 647 Mouse Anti-Human CD16


Product Information

| | |
|------------------|---|
| Material Number: | 557710 |
| Size: | 100 tests |
| Vol. per Test: | 5.0 µl |
| Clone: | 3G8 |
| Isotype: | Mouse IgG1, κ |
| Reactivity: | Human |
| | Reported: Baboon, Rhesus, Cynomolgus |
| Workshop: | IV N409 |
| Storage Buffer: | Aqueous buffered solution containing BSA and ≤0.09% sodium azide. |

Description

Reacts with the 50-65 kDa transmembrane form of IgG Fc Receptor (FcγRIII), a human NK-cell-associated antigen. CD16 is expressed on NK cells as well as macrophages and granulocytes. Reports indicate that CD16 plays a role in signal transduction and NK cell activation. 3G8 antibody blocks the binding of soluble immune complexes to granulocytes.

This antibody is routinely tested by flow cytometric analysis. Other applications were tested at BD Biosciences Pharmingen during antibody development only or reported in the literature.

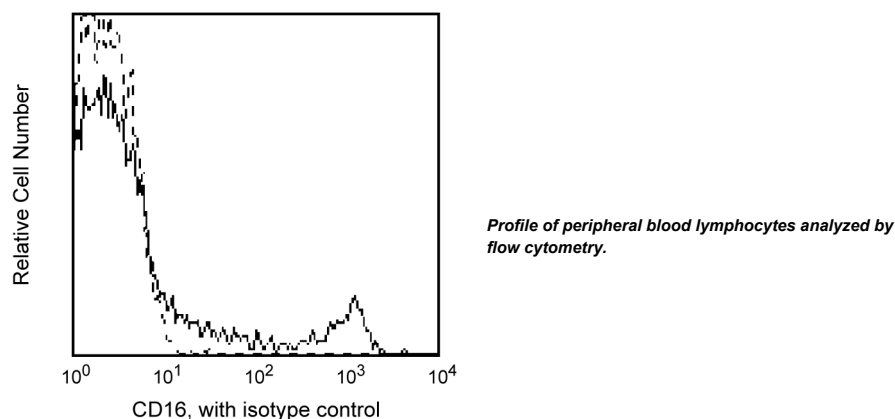

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed. Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

| | |
|----------------|------------------|
| Flow cytometry | Routinely Tested |
|----------------|------------------|

Suggested Companion Products

| Catalog Number | Name | Size | Clone |
|----------------|---|-----------|---------|
| 557714 | Alexa Fluor® 647 Mouse IgG1 κ Isotype Control | 100 tests | MOPC-21 |

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100-µl experimental sample (a test).
2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
3. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

BD Biosciences

bdbiosciences.com

| United States | Canada | Europe | Japan | Asia Pacific | Latin America/Caribbean |
|---------------|--------------|---------------|--------------|--------------|-------------------------|
| 877.232.8995 | 888.259.0187 | 32.53.720.550 | 0120.8555.90 | 65.6861.0633 | 55.11.5185.9995 |

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



4. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
5. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
6. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
7. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
8. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
9. Alexa Fluor is a registered trademark of Molecular Probes, Inc., Eugene, OR.

References

Barclay NA, Brown MH, Birkeland ML, et al, ed. *The Leukocyte Antigen FactsBook*. San Diego, CA: Academic Press; 1997. (Biology)

Fleit HB, Wright SD, Unkeless JC. Human neutrophil Fc gamma receptor distribution and structure. *Proc Natl Acad Sci U S A*. 1982; 79(10):3275-3279. (Biology)

Knapp W, Dorken B, Rieber EP, et al, ed. *Leucocyte Typing IV*. New York: Oxford University Press; 1989. (Clone-specific)

Stroncek DF, Skubitz KM, Plachta LB, et al. Alloimmune neonatal neutropenia due to an antibody to the neutrophil Fc-gamma receptor III with maternal deficiency of CD16 antigen. *Blood*. 1991; 77(7):1572-1580. (Biology)

Wirthmueller U, Kurosaki T, Murakami MS, Ravetch JV.. Signal transduction by Fc gamma RIII (CD16) is mediated through the gamma chain. *J Exp Med*. 1992; 175(5):1381-1390. (Biology)