Technical Data Sheet

Alexa Fluor® 700 Rat IgG1 κ Isotype Control

Product Information

558001 **Material Number:** 0.1 mg Size: 0.2 mg/ml **Concentration:** R3-34 Clone:


Mouse immunoglobulin Immunogen:

Rat IgG1, ĸ Isotype:

Aqueous buffered solution containing protein stabilizer and ≤0.09% sodium Storage Buffer:

Description

The immunoglobulin from the R3-34 hybridoma was identified as a non-reactive clone, following immunization with mouse Ig. The R3-34 immunoglobulin was selected as an isotype control following screening for low background on a variety of mouse and human tissues.

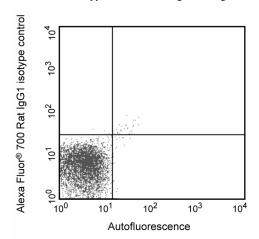

Profile of Rat IgG1 (R3-34) intracellular staining of MiCK-1 cells, Cat. No. 554652 (activated mouse splenocytes) analyzed by flow cytometry

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated to Alexa Fluor® 700 under optimum conditions, and unreacted Alexa Fluor® 700 was removed.

Store undiluted at 4° C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

-	Tippication		
	Intracellular staining (flow cytometry)	Routinely Tested	
	Isotype control	Routinely Tested	

Recommended Assay Procedure:

An isotype control should be used at the same concentration as the antibody of interest.

Suggested Companion Products

Catalog Number	Name	Size	Clone
554715	BD Cytofix/Cytoperm Plus Kit (with BD GolgiStop)	250 tests	(none)
555028	BD Cytofix/Cytoperm Plus Kit (with BD GolgiPlug)	250 tests	(none)

Product Notices

Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

BD Biosciences

www.bdbiosciences.com

Europe 32.53.720.550 0120.8555.90 877.232.8995 888.259.0187 65.6861.0633 55.11.5185.9995 For country-specific contact information, visit www.bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation Conditions: The information disclosed nerein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. @2007 BD



- Alexa Fluor® 700 has an adsorption maximum of ~700nm and a peak fluorescence emission of ~720nm. Before staining cells with this
 reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
- 3. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/pharmingen/colors.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 5. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.

558001 Rev. 4 Page 2 of 2