

## Technical Data Sheet

## Alexa Fluor® 488 Mouse anti-PLCγ1


## Product Information

Material Number:	558566
Size:	50 tests
Vol. per Test:	20 µl
Clone:	10/PLCgamma
Immunogen:	Cow PLCγ1 N-terminal region Peptide
Isotype:	Mouse IgG1
Reactivity:	Confirmed by flow cytometry: Human Confirmed by western blot using purified antibody (Cat. No. 610027 or 610028): Chicken, Dog, Human, Mouse, Rat
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

## Description

The Phospholipase C (PLC) isozymes hydrolyze phosphatidyl inositol biphosphate to inositol triphosphate and diacylglycerol. The former causes release of calcium from the endoplasmic reticulum, while the latter is an activator of Protein Kinase C. Within the PLC family, PLCγ is the only member that contains SH2 and SH3 domains. These domains enable it to interact with receptor tyrosine kinases and become enzymatically activated via phosphorylation. It exists as two isoforms: 1) PLCγ1, which is ubiquitously expressed, and 2) PLCγ2, found primarily in the lymphoid system. PLCγ is essential for growth factor-induced cell motility and mitogenesis. PLCγ1 null mice exhibit retarded embryonic growth and lethality in midgestation. Overexpression of PLCγ is evident in several forms of cancer, and it has been identified as a key mediator of PDGF-dependent cellular transformation. Thus regulation of PLCγ activity by growth factors is involved in cell growth and transformation.

The 10/PLCgamma monoclonal antibody recognizes PLCγ1, regardless of phosphorylation status. It does not cross-react with PLCγ2.

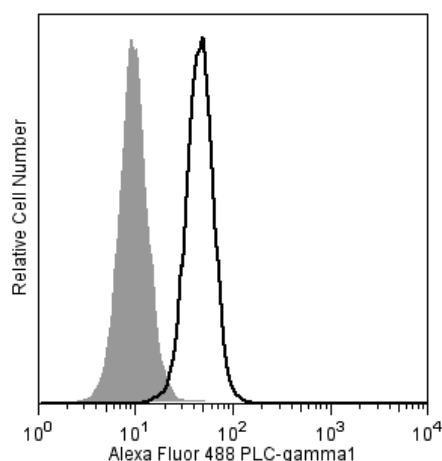

**Analysis of PLCγ1 in human T leukemia cell lines.**  
Jurkat (ATCC TIB-152, open histogram) and J.gamma1 (PLCγ1-deficient Jurkat, ATCC CRL-2678, shaded histogram) were fixed with pre-warmed BD Cytofix™ Fixation buffer (Cat. No. 554655) for 10 minutes at 37°C, permeabilized with BD Phosflow™ Perm Buffer III (Cat. No. 558050) on ice for at least 30 minutes, and then stained with Alexa Fluor® 488 Mouse anti-PLCγ1 (Cat. No. 558566). Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system.

## Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 488 under optimum conditions, and unreacted Alexa Fluor® 488 was removed.

## Application Notes

## Application

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

## Suggested Companion Products

Catalog Number	Name	Size	Clone
558050	Perm Buffer III	125 ml	(none)
554655	Fixation Buffer	100 ml	(none)

## BD Biosciences

[bdbiosciences.com](http://bdbiosciences.com)

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



## Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100- $\mu$ l experimental sample (a test).
2. Please refer to [www.bdbiosciences.com/pharming/en/protocols](http://www.bdbiosciences.com/pharming/en/protocols) for technical protocols.
3. Alexa Fluor® 488 fluorochrome emission is collected at the same instrument settings as for fluorescein isothiocyanate (FITC).
4. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
5. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
7. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
8. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.

## References

Chen P, Murphy-Ullrich JE, Wells A. A role for gelsolin in actuating epidermal growth factor receptor-mediated cell motility. *J Cell Biol.* 1996; 134(3):689-698. (Biology)

Obermeier A, Tinhofer I, Grunicke HH, Ullrich A. Transforming potentials of epidermal growth factor and nerve growth factor receptors inversely correlate with their phospholipase C gamma affinity and signal activation. *EMBO J.* 1996; 15(1):73-82. (Biology)

## BD Biosciences

[bdbiosciences.com](http://bdbiosciences.com)

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

