Technical Data Sheet

APC Mouse Anti-Human CD29

Product Information

Material Number: 559883

Alternate Name: Integrin \$1 Chain

100 tests Size Vol. per Test: 20 ul MAR4 Clone: Isotype: Mouse IgG1, κ Reactivity: QC Testing: Human


Tested in Development: Baboon, Rhesus, Cynomolgus.

Workshop:

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Reacts with the 130 kDa integrin β1 subunit that is expressed as a heterodimeric complex with one of six distinct α subunits, comprising the very late activation antigen (VLA) subfamily of adhesion receptors. The β1 subunit has a broad tissue distribution; it is expressed on lymphocytes, monocytes and weakly on granulocytes, but not on erythrocytes. These receptors are involved in a variety of cell-cell and cell-matrix interactions.

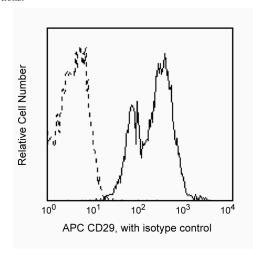

Profile of peripheral blood lymphocytes analyzed by flow cytometry.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number	Name	Size	Clone
555751	APC Mouse IgG1, κ Isotype Control	100 tests	MOPC-21
554681	APC Mouse IgG1 κ Isotype Control	0.1 mg	MOPC-21
554656	Stain Buffer (FBS)	500 ml	(none)

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

BD Biosciences

bdbiosciences.com

Asia Pacific 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



559883 Rev. 6

- 4. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 6. This APC-conjugated reagent can be used in any flow cytometer equipped with a dye, HeNe, or red diode laser.

References

Hemler ME. VLA proteins in the integrin family: structures, functions, and their role on leukocytes. *Annu Rev Immunol.* 1990; 8:365-400. (Biology) Hynes RO. Integrins: versatility, modulation, and signaling in cell adhesion. *Cell.* 1992; 69(1):11-25. (Biology)

Kishimoto T, von dem Borne AEG, Goyert SM,et al., ed. *Leucocyte Typing VI: White Cell Differentiation Antigens*. London: Garland Publishing; 1997. (Clone-specific)

559883 Rev. 6 Page 2 of 2