

Technical Data Sheet

Alexa Fluor® 488 Mouse anti-PKCα (pT497)


Product Information

Material Number:	560139
Size:	50 Tests
Vol. per Test:	20 µl
Clone:	K14-984
Immunogen:	Phosphorylated Human PKCα Peptide
Isotype:	Mouse (BALB/c) IgG1, κ
Reactivity:	QC Testing: Human Tested in Development: Mouse
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The Protein Kinase C (PKC) family of serine/threonine protein kinases is involved in a number of processes such as growth, differentiation, and cytokine secretion. Three categories exist, conventional PKC (cPKC), novel PKC (nPKC), and atypical PKC (aPKC). All have C-terminal kinase domains, which are closely related to those of protein kinases A and B (Akt), and variable N-terminal regulatory domains that give them different modes of activation. For example, cPKC (α, βI, βII, and γ isoforms) are calcium-activated, phospholipid-dependent serine/threonine-specific enzymes that can also be activated by phorbol esters. However, nPKC (δ, ε, η, and θ isoforms) and aPKC (ζ, ι, and λ isoforms) are Ca²⁺-independent. aPKC are unique in that their activity is independent of diacylglycerols and phorbol esters. Phosphorylation at three conserved sites in the kinase domain is required for catalytic activity. Specifically, the threonine 497 (T497) of PKCα is in the activation loop of the kinase domain and is phosphorylated by the constitutively active phosphoinositide-dependent kinase-1 (PDK-1).

The K14-984 monoclonal antibody recognizes the phosphorylated T497 in the kinase domain of human PKCα.

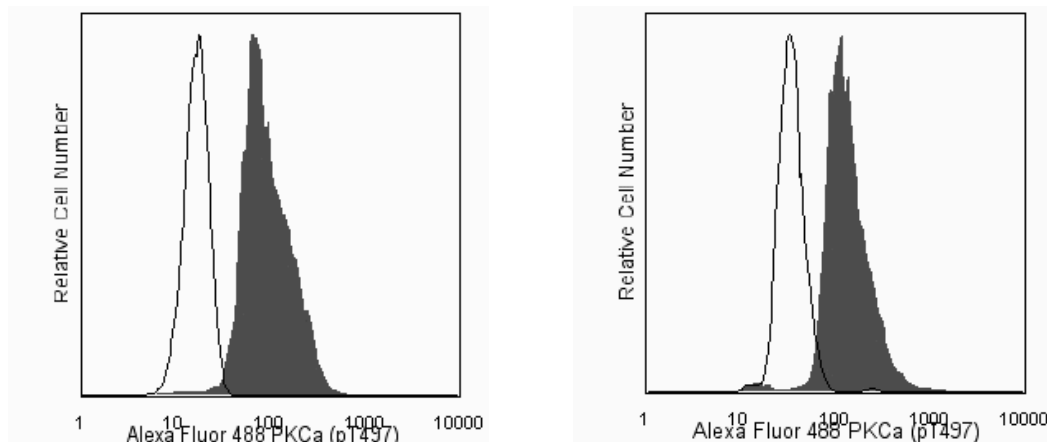

Analysis of PKCα (pT497) in human peripheral blood lymphocytes and vascular endothelium.

LEFT PANEL: Human peripheral blood mononuclear cells were either treated with 50 nM calyculin A for 30 minutes at 37°C (shaded histogram) or untreated (open histogram). For data analysis, lymphocytes were selected by their scatter profile.

RIGHT PANEL: After serum starvation overnight, EA.hy926 cells (see reference Edgell, McDonald, Graham, 1983; ATCC CRL-2922) were detached with trypsin, washed, resuspended in serum-free DMEM, and either treated with 50 nM calyculin A (shaded histogram) for 30 minutes at 37°C or untreated (open histogram).

The cells were fixed (BD Cytofix™ buffer, Cat. No. 554655) for 10 minutes at 37°C, permeabilized with BD Phosflow™ Perm Buffer III (Cat. No. 558050) on ice for 30 minutes, and then stained with Alexa Fluor® 488 Mouse anti-PKCα (pT497). The data demonstrates that the level of phosphorylation of PKCα increases when phosphatase activity is inhibited by the treatment. Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system. BD Phosflow™ Fix Buffer I (Cat. No. 557870) may be used for cell fixation.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 488 under optimum conditions, and unreacted Alexa Fluor® 488 was removed.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



The purified or conjugated mAb was characterized by flow cytometry (Flow) and western blot (WB) using these model systems:

Method	Species	Cells	Treatment	Fixation	Perm buffer	Result
Flow	Human	EA-hy 926	Calyculin A	Cytofix or Fix I	Perm III	Up-regulation
	Human	PBMC	Calyculin A	Cytofix or Fix I	Perm III	Up-regulation
	Human	Jurkat	Phospho peptide	Cytofix or Fix I	Perm III	Blocking
	Human	Jurkat	Irrelevant peptide	Cytofix or Fix I	Perm III	No blocking
	Human	Jurkat	Lambda phosphatase	Cytofix or Fix I	Perm III	Signal reduction
	Human	Jurkat	Kinase inhibitors	Cytofix or Fix I	Perm III	No change
WB	Human	EA-hy 926	Calyculin A			Up-regulation
	Human	PBMC	Calyculin A			Up-regulation
	Human	Jurkat	Lambda phosphatase			Signal reduction
	Mouse	Splenocytes	Calyculin A			Up-regulation

Application Notes

Application

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

Suggested Companion Products

Catalog Number	Name	Size	Clone
554655	Fixation Buffer	100 mL	(none)
558050	Perm Buffer III	125 mL	(none)
557870	Fix Buffer I	250 mL	(none)
554656	Stain Buffer (FBS)	500 mL	(none)

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. Please refer to www.bdbiosciences.com/pharming/en/protocols for technical protocols.
3. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
4. Alexa Fluor® 488 fluorochrome emission is collected at the same instrument settings as for fluorescein isothiocyanate (FITC).
5. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
7. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
8. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
9. Species testing during development may have been performed with a different format of the same clone. Selected applications have been tested for cross-reactivity.

References

Edgell C-JS, McDonald CC, Graham JB. Permanent cell line expressing human factor VIII-related antigen established by hybridization. *Proc Natl Acad Sci U S A*. 1983; 80:3734-3737. (Methodology)

Newton AC. Regulation of the ABC kinases by phosphorylation: protein kinase C as a paradigm. *Biochem J*. 2003; 370:361-371. (Biology)

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 800.268.5430 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

