

Technical Data Sheet


Alexa Fluor® 555 Mouse anti-Oct3/4

Product Information

Material Number:	560306
Alternate Name:	Oct3, OTF3, Oct4, OTF4, POU5F1
Size:	100 tests
Vol. per Test:	5 µl
Clone:	40/Oct-3
Immunogen:	Mouse Oct3 aa. 252-372 Recombinant Protein
Isotype:	Mouse IgG1, κ
Reactivity:	QC Testing: Human Tested in development using purified antibody: Mouse
Storage Buffer:	Aqueous buffered solution containing BSA, protein stabilizer, and ≤0.09% sodium azide.

Description

Development of a multicellular organism from a single fertilized cell is regulated by the coordinated activity of DNA transcription factors. Oct3/4, a member of the POU family of transcription factors, functions in pluripotent cells of early embryonic stem cell (ES) lines and embryonal carcinomas (EC). Other members of the POU family include Oct1, Oct2, Pit-1, and unc-86. The POU domain, a 150-amino acid region that determines binding specificity, is conserved among these proteins and consists of 3 subdomains: POU-specific A and B subdomains and a homeobox-like subdomain. Oct3/4 is expressed in undifferentiated cells, but is lost as cells are induced to differentiate. Oct3/4 is not expressed in adult tissues. The interaction of Oct3/4 with SOX2, another embryonic transcription factor, produces an active complex that regulates expression of genes such as Nanog, UTF1, and FGF4. Although Oct3/4 is specifically phosphorylated on serine residues, this modification is not required for DNA binding, but may affect its transactivation potential. Thus, Oct3/4 is a transcription factor that plays an important role in determining early steps of embryogenesis and differentiation.

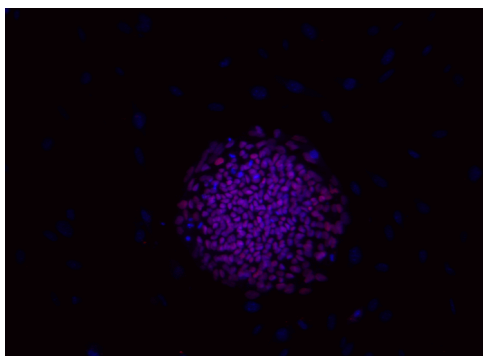

Immunofluorescent staining of human ES cell line. H9 cells (WiCell, Madison, WI) were grown on mouse epithelial feeder cells in a BD Falcon™ 96-well Imaging Plate (Cat. No. 353219). After overnight incubation, the cells were fixed, permeabilized with BD Perm/Wash™ Buffer (Cat. No. 554723), stained with Alexa Fluor® 555 Mouse anti-Oct 3/4 (pseudo colored red), and counter stained with Hoechst 33342 (pseudo colored blue) according to the Recommended Assay Procedure. The images were captured on a BD Pathway™ 435 High-Content Bioimager System with a 20x objective and merged using BD AttoVision™ software. This antibody also stained human EC NCCIT (ATCC, Cat. No. CRL-2073), mouse EC F9 (ATCC, Cat. No. CRL-1720), and mouse ES-E14TG2a (ATCC, Cat. No. CRL-1821) cells. It also works with Triton™ X-100 and methanol permeabilization (see Recommended Assay Procedure).

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated to Alexa Fluor® 555 under optimum conditions, and unreacted Alexa Fluor® 555 was removed. Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

Bioimaging	Routinely Tested
------------	------------------

Recommended Assay Procedure:

For more information, please refer to: Bioimaging: http://www.bdbiosciences.com/pharmingen/protocols/Bioimaging_Certified.shtml

Recommended Protocol for Bioimaging

- Seed the cells in appropriate culture medium at an appropriate cell density in a BD Falcon™ 96-well Imaging Plate (Cat. No. 353219), and culture overnight to 48 hours.
- Remove the culture medium from the wells, and wash (one to two times) with 100 µl of 1× PBS.
- Fix the cells by adding 100 µl of fresh 3.7% Formaldehyde in PBS or BD Cytotfix™ fixation buffer (Cat. No. 554655) to each well and incubating for 10 minutes at room temperature (RT).

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



4. Remove the fixative from the wells, and wash the wells (one to two times) with 100 µl of 1× PBS.
5. Permeabilize the cells using either cold methanol (a), Triton™ X-100 (b), or Saponin (c):
 - a. Add 100 µl of -20°C 90% methanol or -20°C BD™ Phosflow Perm Buffer III (Cat. No. 558050) to each well and incubate for 5 minutes at RT.
 - b. Add 100 µl of 0.1% Triton™ X-100 to each well and incubate for 5 minutes at RT.
 - c. Add 100 µl of 1× Perm/Wash buffer (Cat. No. 554723) to each well and incubate for 15 to 30 minutes at RT. Continue to use 1× Perm/Wash buffer for all subsequent wash and dilutions steps.
6. Remove the permeabilization buffer from the wells, and wash one to two times with 100 µl of appropriate buffer (either 1× PBS or 1× Perm/Wash buffer, see step 5.c.).
7. Optional blocking step: Remove the wash buffers, and block the cells by adding 100 µl of blocking buffer BD Pharmingen™ Stain Buffer (FBS) (Cat. No. 554656) or 3% FBS in appropriate dilution buffer to each well and incubating for 15 to 30 minutes at RT.
8. Dilute the antibody to its optimal working concentration in appropriate dilution buffer. Titrate purified (unconjugated) antibodies and second-step reagents to determine the optimal concentration. If using a Bioimaging Certified antibody conjugate, dilute it 1:10.
9. Add 50 µl of diluted antibody per well and incubate for 60 minutes at RT. Incubate in the dark if using fluorescently labeled antibodies.
10. Remove the antibody, and wash the wells three times with 100 µl of wash buffer. An optional detergent wash (100 µl of 0.05% Tween in 1× PBS) can be included prior to the regular wash steps.
11. If the antibody being used is fluorescently labeled, then move to step 12. Otherwise, if using a purified unlabeled antibody, repeat steps 8 to 10 with a fluorescently labeled second-step reagent to detect the purified antibody.
12. After the final wash, counter-stain the nuclei by adding 100 µl of a 2 g/ml solution of Hoechst 33342 (eg, Sigma-Aldrich Cat. No. B2261) in 1× PBS to each well at least 15 minutes before imaging.
13. View and analyze the cells on an appropriate imaging instrument. Recommended filters for the BD Pathway™ instruments are:

Instrument	Excitation	Emission	Dichroic
BD Pathway 855	548/20	570LP	Fura/FITC
BD Pathway 435	543/22	593/40	FF562

Suggested Companion Products

Catalog Number	Name	Size	Clone
353219	BD Falcon™ 96-well Imaging Plate	NA	(none)
554723	Perm/Wash Buffer	100 ml	(none)
554655	Fixation Buffer	100 ml	(none)
554656	Stain Buffer (FBS)	500 ml	(none)
558050	Perm Buffer III	125 ml	(none)
611202	Purified Mouse Anti-Oct3/4	50 µg	40/Oct-3
611203	Purified Mouse Anti-Oct-3/4	150 µg	40/Oct-3

Product Notices

1. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
2. This reagent has been pre-diluted for use at the recommended Volume per Test when following the Recommended Assay Procedure. A Test is typically ~10,000 cells cultured in a well of a 96-well imaging plate.
3. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
4. Triton is a trademark of the Dow Chemical Company.
5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
6. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
7. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.

References

Nishimoto M, Fukushima A, Okuda A, Muramatsu M. The gene for the embryonic stem cell coactivator UTF1 carries a regulatory element which selectively interacts with a complex composed of Oct-3/4 and Sox-2. *Mol Cell Biol.* 1999; 19(8):5453-5465. (Biology)

Okamoto K, Okazawa H, Okuda A, Sakai M, Muramatsu M, Hamada H. A novel octamer binding transcription factor is differentially expressed in mouse embryonic cells. *Cell.* 1990; 60(3):461-472. (Biology)

Pan G, Thomson JA. Nanog and transcriptional networks in embryonic stem cell pluripotency. *Cell Res.* 2007; 17:42-49. (Biology)

Rosfjord E, Scholtz B, Lewis R, Rizzino A. Phosphorylation and DNA binding of the octamer binding transcription factor Oct-3. *Biochem Biophys Res Commun.* 1995; 212(3):847-853. (Biology)

Vigano MA, Staudt LM. Transcriptional activation by Oct-3: evidence for a specific role of the POU-specific domain in mediating functional interaction with Oct-1. *Nucleic Acids Res.* 1996; 24(11):2112-2118. (Biology)

Yuan H, Corbi N, Basilico C, Dailey L. Developmental-specific activity of the FGF-4 enhancer requires the synergistic action of Sox2 and Oct-3. *Genes Dev.* 1995; 9(21):2635-2645. (Biology)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

