Technical Data Sheet

V450 Mouse Anti-Human CD3

Product Information

Reactivity: Human

QC Testing: Baboon or Rhesus or Cynomolgus

Workshop: NA


Storage Buffer: Aqueous buffered solution containing protein stabilizer and ≤0.09% sodium

azide.

Description

Clone SP34-2 is a mouse IgG1 isotype monoclonal antibody, descendant of SP34 (mouse IgG3), with the same specificity and reactivity pattern as the parent clone. It cross-reacts with a major subset of peripheral blood lymphocytes, but not monocytes or granulocytes, of baboon, and rhesus, cynomolgus, and pigtail macaque monkeys. The distribution on lymphocytes is similar to that observed with normal human donor lymphocytes with the majority of CD3-positive cells being negative when dual stained with antibodies to B or NK cells markers. SP34-2 is also capable of inducing cell proliferation on both human and non-human primate PBMC.

The antibody is conjugated to BD HorizonTM V450, which has been developed for use in multicolor flow cytometry experiments and is available exclusively from BD Biosciences. It is excited by the Violet laser Ex max of 406 nm and has an Em Max at **450** nm. Conjugates with BD HorizonTM V450 can be used in place of Pacific BlueTM conjugates.

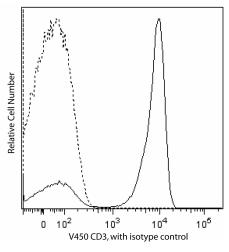

Analysis of CD3 on peripheral blood lymphocytes of rhesus macaque (Macaca mulatta) monkey. Peripheral blood lymphocytes were stained with BD Horizon™ V450 Mouse Anti-Human CD3 (solid line) and compared to whole blood stained with BD Horizon™ V450 Mouse IgG1, κ Isotype Control (clone MOPC-21, Cat. No. 560373, dotted line). Lymphocytes were selected by light scatter profile. Flow cytometry was performed on a BD LSR™ If flow cytometry system.

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD HorizonTM V450 under optimum conditions, and unreacted BD HorizonTM V450 was removed. Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

•	-phieuron				
[]	Flow cytometry	Routinely Tested			

BD Biosciences

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 888.259.0187
 32.53.720.550
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale. BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



Suggested Companion Products

Catalog Number	Name	Size	Clone
560373	V450 Mouse IgG1, κ Isotype Control	0.1 mg	MOPC-21

Product Notices

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁶ cells in a 100-μl experimental sample (a test).
- 2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 3. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 4. BD HorizonTM V450 has a maximum absorption of 406 nm and maximum emission of 450 nm. Before staining with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
- Pacific Blue™ is a trademark of Molecular Probes, Inc., Eugene, OR.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

References

Schlossman S, Boumell L, et al, ed. Leucocyte Typing V. New York: Oxford University Press; 1995.(Clone-specific: Induction)

Bleavins MR, Brott DA, Alvey JD, de la Iglesia FA. Flow cytometric characterization of lymphocyte subpopulations in the cynomolgus monkey (Macaca fascicularis). Vet Immunol Immunopathol. 1993; 37(1):1-13.(Biology)

Conrad ML, Davis WC, Koop BF. TCR and CD3 antibody cross-reactivity in 44 species. Cytometry A. 2007; 71(11):925-933.(Biology)

Engel P, Wagner N, Zhou L, et al. CD86 Workshop Report. In: Schlossman SF, Boumsell L, Gilks W, et al, ed. Leukocyte Typing V: White Cell Differentiation Antigens. New York: Oxford University Press; 1995.(Biology)

Jacobsen CN, Aasted B, Broe MK, Petersen JL. Reactivities of 20 anti-human monoclonal antibodies with leucocytes from ten different animal species. *Vet Immunol Immunopathol.* 1993; 39(4):461-466.(Biology)

Wilson AD, Shooshtari M, Finerty S, Watkins P, Morgan AJ. Selection of monoclonal antibodies for the identification of lymphocyte surface antigens in the New World primate Saguinus oedipus oedipus (cotton top tamarin). *J Immunol Methods*. 1995; 178(2):195-200.(Biology)

560352 Rev. 1 Page 2 of 2