


Technical Data Sheet

Alexa Fluor® 700 Mouse Anti-Human CD45**Product Information**

Material Number:	560566
Alternate Name:	Leukocyte Common Antigen
Size:	50 tests
Vol. per Test:	5 µl
Clone:	HI30
Isotype:	Mouse IgG1, κ
Reactivity:	QC Testing: Human
Workshop:	IV N816
Storage Buffer:	Aqueous buffered solution containing protein stabilizer and ≤0.09% sodium azide.

Description

Reacts with the 180, 190, 205, 220 kD isoforms of the leukocyte common antigen (LCA) present on all human leukocytes including lymphocytes, monocytes, granulocytes, eosinophils, and thymocytes but absent from circulating erythrocytes, platelets, or mature erythroid cells of bone marrow and non-hemopoietic tissues.

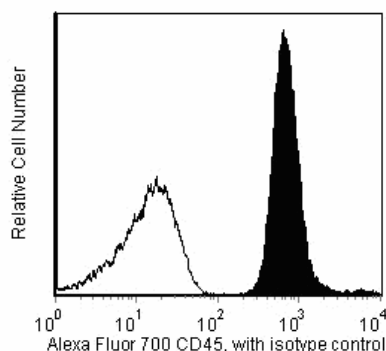

Flow cytometric analysis of CD45 on human lysed whole blood. Human lysed whole blood was stained with the Alexa Fluor® 700 Mouse Anti-Human CD45 antibody (shaded) or with a Alexa Fluor® 700 Mouse IgG1,κ isotype control (unshaded). Histograms were derived from gated events based on light scattering characteristics for lymphocytes. Flow cytometry was performed on a BD™ LSR II flow cytometry system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 700 under optimum conditions, and unreacted Alexa Fluor® 700 was removed.

Application Notes**Application**

Flow cytometry	Routinely Tested
----------------	------------------

Suggested Companion Products

Catalog Number	Name	Size	Clone
557882	Alexa Fluor® 700 Mouse IgG1, κ Isotype Control	0.1 mg	MOPC-21

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100-µl experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Alexa Fluor® 700 has an adsorption maximum of ~700nm and a peak fluorescence emission of ~720nm. Before staining cells with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
4. Alexa Fluor is a registered trademark of Molecular Probes, Inc., Eugene, OR.
5. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
7. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
8. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.

References

- Biffen M, McMichael-Phillips D, Larson T, Venkitaraman A, Alexander D. The CD45 tyrosine phosphatase regulates specific pools of antigen receptor-associated p59fyn and CD4-associated p56lck tyrosine in human T-cells. *EMBO J.* 1994; 13(8):1920-1929. (Biology)
- Hermiston ML, Xu Z, Weiss A. CD45: a critical regulator of signaling thresholds in immune cells. *Annu Rev Immunol.* 2003; 21:107-137. (Biology)
- Knapp W, Dorken B, Rieber EP, et al, ed. *Leucocyte Typing IV*. New York: Oxford University Press; 1989. (Biology)
- Loken MR, Brosnan JM, Bach BA, Ault KA. Establishing optimal lymphocyte gates for immunophenotyping by flow cytometry. *Cytometry.* 1990; 11(4):453-459. (Methodology: Flow cytometry)
- Phipps DJ, Yousefi S, Branch DR. Increased enzymatic activity of the T-cell antigen receptor-associated fyn protein tyrosine kinase in asymptomatic patients infected with the human immunodeficiency virus. *Blood.* 1997; 90(9):3603-3612. (Biology)
- Sieh M, Bolen JB, Weiss A. CD45 specifically modulates binding of Lck to a phosphopeptide encompassing the negative regulatory tyrosine of Lck. *EMBO J.* 1994; 12(1):315-321. (Biology)
- Terstappen LW, Levin J. Bone marrow cell differential counts obtained by multidimensional flow cytometry. *Blood Cells.* 1992; 18(2):311-330. (Biology)
- Trowbridge IS, Thomas ML. CD45: an emerging role as a protein tyrosine phosphatase required for lymphocyte activation and development. *Annu Rev Immunol.* 1994; 12:85-116. (Biology)
- Wang Y, Liang L, Esselman WJ. Regulation of the calcium/NF-AT T cell activation pathway by the D2 domain of CD45. *J Immunol.* 2000; 164(5):2557-2564. (Biology)
- Wu L, Fu J, Shen SH. SKAP55 coupled with CD45 positively regulates T-cell receptor-mediated gene transcription. *Mol Biol Cell.* 2002; 22(8):2673-2686. (Biology)