Technical Data Sheet

Alexa Fluor® 700 Mouse Anti-Human CD20

Product Information

Material Number: 560631 50 tests Size: 5 µl Vol. per Test: 2H7 Clone:


Mouse IgG2b, κ Isotype: QC Testing: Human Reactivity:


IV B201 Workshop:

Storage Buffer: Aqueous buffered solution containing protein stabilizer and ≤0.09% sodium

Description

Reacts with a 33-37 kDa unglycosylated phosphoprotein which crosses the cell membrane. CD20 is expressed on pre-B-cells, resting and activated B cells and follicular dendritic cells but not plasma cells. Low level CD20 expression is observed on a small subset of normal circulating T lymphocytes. The CD20 molecule is involved in regulation of B-cell activation.

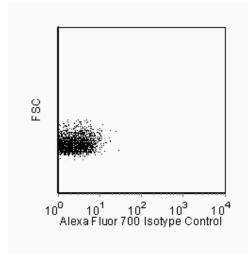

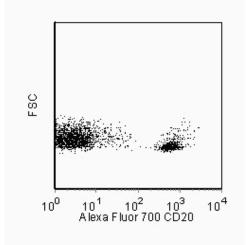

Flow cytometric analysis of CD20 on human lysed whole blood. Human lysed whole blood was stained either with a Alexa Fluor® 700 Mouse IgG2b, κ isotype control (left panel) or with the Alexa Fluor® 700 Mouse Anti-Human CD20 antibody (right panel). Dot plots were derived from gated events based on light scattering characteristics for lymphocytes. Flow cytometry was performed on a BD™ LSR II flow cytometry system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 700 under optimum conditions, and unreacted Alexa Fluor® 700 was removed.

Application Notes

Application

_	r		
	Flow cytometry	Routinely Tested	

Suggested Companion Products

Catalog Number	Name	Size	Clone
560543	Alexa Fluor® 700 Mouse IgG2b, κ Isotype Control	0.1 mg	27-35
555899	Lysing Buffer	100 ml	(none)

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁶ cells in a 100-µl experimental sample (a test).
- An isotype control should be used at the same concentration as the antibody of interest.
- Alexa Fluor® 700 has an adsorption maximum of ~700nm and a peak fluorescence emission of ~720nm. Before staining cells with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
- Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.

BD Biosciences

bdbiosciences.com **United States** Canada Asia Pacific Latin America/Caribbean Europe 877.232.8995 888.259.0187 32.53.720.550 0120.8555.90 65.6861.0633 55.11.5185.9995 For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



- 5. The Alexa Fluor®, Pacific BlueTM, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific BlueTM dye, and Cascade Blue® dye are covered by pending and issued patents.
- 6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 7. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 8. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

Hultin LE, Hausner MA, Hultin PM, Giorgi JV. CD20 (pan-B cell) antigen is expressed at a low level on a subpopulation of human T lymphocytes. Cytometry. . (Biology)

Knapp W, Dorken B, Rieber EP, et al, ed. *Leucocyte Typing IV*. New York: Oxford University Press; 1989. (Biology) Schlossman S, Boumell L, et al, ed. *Leucocyte Typing V*. New York: Oxford University Press; 1995. (Biology)

560631 Rev. 1 Page 2 of 2