## Technical Data Sheet

# APC Mouse Anti-Human CD2

#### **Product Information**


**Material Number:** 560642 50 tests Size: 5 µl Vol. per Test: RPA-2.10 Clone: Mouse IgG1, κ Isotype: QC Testing: Human Reactivity:

IV T085 Workshop:

Aqueous buffered solution containing BSA and ≤0.09% sodium azide. Storage Buffer:

#### Description

CD2 reacts with the 50 kDa single-chain transmembrane glycoprotein, also known as LFA-2 or the receptor for sheep erythrocytes. CD2 belongs to the immunoglobulin superfamily of proteins along with its ligand LFA-3, CD58. It is present on about 80-90% of human peripheral blood lymphocytes, greater than 95% of thymocytes, all T lymphocytes that form E-rosettes and a subset of NK cells. CD2 plays a role in T-cell signaling and in lymphocyte adhesion.

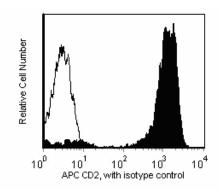

Flow cytometric analysis of CD2 on human lysed whole blood. Human lysed whole blood was stained with the APC Mouse Anti-Human CD2 antibody (shaded) or with a APC Mouse IgG1, κ isotype control (unshaded). Histograms were derived from gated events based on light scattering characteristics for lymphocytes. Flow cytometry was performed on a BD™ LSR II flow cytometry system.

## **Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

#### **Application Notes**

Flow cytometry Routinely Tested

## **Suggested Companion Products**

Catalog Number	Name	Size	Clone
555751	APC Mouse IgG1, κ Isotype Control	100 tests	MOPC-21
555899	Lysing Buffer	100 ml	(none)

#### **Product Notices**

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100- $\mu$ l experimental sample (a test).
- 2. An isotype control should be used at the same concentration as the antibody of interest.
- 3 This APC-conjugated reagent can be used in any flow cytometer equipped with a dye, HeNe, or red diode laser.
- Source of all serum proteins is from USDA inspected abattoirs located in the United States. 4.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

## **BD Biosciences**

bdbiosciences.com

**United States** Asia Pacific Latin America/Caribbean Europe 877.232.8995 888.259.0187 32.53.720.550 0120.8555.90 65.6861.0633 55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how\_to\_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



560642 Rev. 1

References

Hahn WC, Burakoff SJ, Bierer BE. Signal transduction pathways involved in T cell receptor-induced regulation of CD2 avidity for CD58. *J Immunol.* 1993; 150(7):2607-2619. (Biology)

Knapp W, Dorken B, Rieber EP, et al, ed. *Leucocyte Typing IV*. New York: Oxford University Press; 1989. (Biology)

Sopper S, Stahl-Hennig C, Demuth M, Johnston IC, Dorries R, ter Meulen V. Lymphocyte subsets and expression of differentiation markers in blood and lymphoid organs of rhesus monkeys. *Cytometry.* 1997; 29(4):351-362. (Biology)

560642 Rev. 1 Page 2 of 2