
Technical Data Sheet APC-H7 Mouse Anti-Human CD20

Product Information

Material Number:	560734
Size:	50 tests
Vol. per Test:	5 μl
Clone:	2H7
Isotype:	Mouse IgG2b, κ
Reactivity:	QC Testing: Human
Workshop:	IV B201
Storage Buffer:	Aqueous buffered solution containing BSA, protein stabilizer, and $\leq 0.09\%$
	sodium azide.

Description

Reacts with a 33-37 kDa unglycosylated phosphoprotein which crosses the cell membrane. CD20 is expressed on pre-B-cells, resting and activated B cells and follicular dendritic cells but not plasma cells. Low level CD20 expression is observed on a small subset of normal circulating T lymphocytes. The CD20 molecule is involved in regulation of B-cell activation.

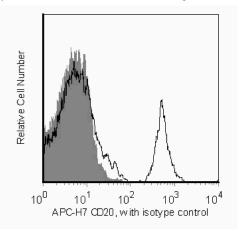

Flow cytometric analysis of CD20 on human lysed whole blood. Human lysed whole blood was stained with the APC-H7 Mouse Anti-Human CD20 antibody (unshaded) or with a APC-H7 Mouse IgG2b, κ isotype control (shaded). Histograms were derived from gated events based on light scattering characteristics for lymphocytes. Flow cytometry was performed on a BD™ LSR II flow cytometry system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with APC-H7 under optimum conditions, and unconjugated antibody and APC-H7 were removed.

Application Notes

Application								
Flow cytometry	Routinely Tested							
Suggested Companion Products								
Catalog Number	Name		Size	Clone				
560183	APC-H7 Mouse IgG2b, κ Isotype Control		0.1 mg	27-35				
555899	Lysing Buffer		100 ml	(none)				

Product Notices

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^{6} cells in a 100-µl experimental sample (a test).
- 2. An isotype control should be used at the same concentration as the antibody of interest.
- 3. Please observe the following precautions: Absorption of visible light can significantly alter the energy transfer occurring in any tandem fluorochrome conjugate; therefore, we recommend that special precautions be taken (such as wrapping vials, tubes, or racks in aluminum foil) to prevent exposure of conjugated reagents, including cells stained with those reagents, to room illumination.

BD Biosciences

bdbiosciences.	com				
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbear
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995
For country-spe	ecific contact int	formation, visit	bdbiosciences.co	m/how_to_orde	r/
of any patents. BL use of our produc product or as a co written authoriza For Research Use	D Biosciences will n its. Purchase does n imponent of anoth tion of Becton Dick Only. Not for use ir	ot be held responsi ot include or carry er product. Any us kinson and Compan diagnostic or thera	ble for patent infrin any right to resell or e of this product oth y is strictly prohibite apeutic procedures.	gement or other vio r transfer this produc ner than the permitte d.	e the above product in violation lations that may occur with the ct either as a stand-alone ed use without the express 108 BD





4. BD APC-H7 is a tandem conjugate and an analog of APC-Cy7 with the same spectral properties. It has decreased intensity but it is engineered for greater stability and less spillover in the APC channel and consequently offers better performance than APC-Cy7. It has an absorption maximum of approximately 650 nm. When excited by light from a red laser, the APC fluorochrome can transfer energy to the cyanine dye, which then emits at a longer wavelength. The resulting fluorescent emission maximum is approximately 767 nm. BD recommends that a 750-nm longpass filter be used along with a red-sensitive detector such as the Hamamatsu R3896 PMT. As with APC-Cy7 special filters are required when using APC-H7 in conjunction with APC.

Note: Although our APC-H7 products demonstrate higher lot-to lot consistency than other APC tandem conjugate products, and every effort is made to minimize the lot-to-lot variation in residual emission from APC, it is strongly recommended that every lot be tested for differences in the amount of compensation required and that individual compensation controls are run for each APC-H7 conjugate. Note: Cy is a trademark of Amersham Biosciences Limited.

- 5. Although BD APC-H7 is engineered to minimize spillover to the APC channel and is more stable and less affected by light, temperature, and formaldehyde-based fixatives, compared to other APC-cyanine tandem dyes, it is still good practice to minimize as much as possible, any light, temperature and fixative exposure when working with all fluorescent conjugates.
- 6. Cy is a trademark of Amersham Biosciences Limited.
- 7. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 8. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 9. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 10. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

Hultin LE, Hausner MA, Hultin PM, Giorgi JV. CD20 (pan-B cell) antigen is expressed at a low level on a subpopulation of human T lymphocytes. Cytometry. . (Biology)

Knapp W, Dorken B, Rieber EP, et al, ed. *Leucocyte Typing IV*. New York: Oxford University Press; 1989. (Biology) Schlossman S, Boumell L, et al, ed. *Leucocyte Typing V*. New York: Oxford University Press; 1995. (Biology)