# **Technical Data Sheet**

# V450 Mouse Anti-Human CD4

#### **Product Information**

560811 **Material Number:** 50 tests Size: 5 µl Vol. per Test: L200 Clone: Isotype: Mouse IgG1, κ


Reactivity: Human

QC Testing: Rhesus or Baboon or Cynomolgus

Storage Buffer: Aqueous buffered solution containing protein stabilizer and ≤0.09% sodium

### Description

The L200 monoclonal antibody specifically binds to the human form of the 56 kDa transmembrane glycoprotein, CD4, present on the T-helper/inducer subset of normal human donor peripheral blood lymphocytes. The L200 antibody also crossreacts with a subset of CD3-positive peripheral blood lymphocytes, but not monocytes, of both Rhesus and Cynomolgus Macaque monkeys. Crossreactivity on both lymphocytes and monocytes (weak) from Baboons is also observed. The distribution on lymphocytes is similar for both human and monkey cells, with the majority of CD4-positive lymphocytes being CD8-negative and lacking reactivity with antibodies to B- or NK-cell markers

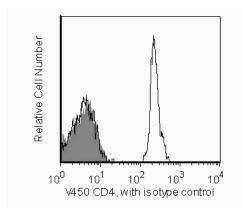

Flow cytometric analysis of CD4 on Rhesus macaque lysed whole blood. Rhesus macaque lysed whole blood was stained with the BD Horizon™ V450 Mouse Anti-Human CD4 antibody (unshaded) or with a BD Horizon™ V450 Mouse IgG1, κ isotype control (shaded). Histograms were derived from gated events based on light scattering characteristics for lymphocytes. Flow cytometry was performed on a BD™ LSR II flow cytometry system.

## **Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon<sup>TM</sup> V450 under optimum conditions, and unreacted BD Horizon<sup>TM</sup> V450 was removed.

## **Application Notes**

### Application

Flow cytometry Routinely Tested

# Suggested Companion Products

Catalog Number	Name	Size	Clone
560373	V450 Mouse IgG1, κ Isotype Control	0.1 mg	MOPC-21
555899	Lysing Buffer	100 ml	(none)

### **Product Notices**

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10<sup>6</sup> cells in a 100-μl experimental
- An isotype control should be used at the same concentration as the antibody of interest.
- BD Horizon<sup>TM</sup> V450 has a maximum absorption of 406 nm and maximum emission of 450 nm. Before staining with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.

## **BD Biosciences**

bdbiosciences.com

**United States** Asia Pacific 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how\_to\_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.
For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding
  to avoid accumulation of potentially explosive deposits in plumbing.
- 5. Pacific Blue<sup>TM</sup> is a trademark of Molecular Probes, Inc., Eugene, OR.
- 6. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 7. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

#### References

Bleavins MR, Brott DA, Alvey JD, de la Iglesia FA. Flow cytometric characterization of lymphocyte subpopulations in the cynomolgus monkey (Macaca fascicularis). Vet Immunol Immunopathol. 1993; 37(1):1-13. (Biology)

Giorgi JV, Hultin LE, Desrosiers RC. The immunopathogenesis of retroviral diseases: no immunophenotypic alterations in T, B, and NK cell subsets in SIVmac239-challenged rhesus macaques protected by SIV delta nef vaccination. *J Med Primatol.* 1996; 25(3):186-191. (Biology)

Indzhiia LV, Yakovleva LA, Overbaugh J, et al. Baboon T cell lymphomas expressing the B cell-associated surface proteins CD40 and Bgp95. *J Clin Invest.* 1992; 12(3):225-236. (Biology)

Jacobsen CN, Aasted B, Broe MK, Petersen JL. Reactivities of 20 anti-human monoclonal antibodies with leucocytes from ten different animal species. *Vet Immunol Immunopathol.* 1993; 39(4):461-466. (Biology)

Powell JD, McClure HM, Anderson D, Fultz PN, Sell KW, Ahmed-Ansari A. Phenotypic and functional differences in NK and LAK cells in the peripheral blood of sooty mangabeys and rhesus macaques. *Cell Immunol.* 1989; 124(1):107-118. (Biology)

Savary CA, Lotzova E, Jackson HJ, Jardine JH, Ang KK. Analysis of interleukin-2-activated killer cells of rhesus monkeys: striking resemblance to the human system. *J Leukoc Biol.* 1993; 54(4):307-313. (Biology)

Schlossman S, Boumell L, et al, ed. Leucocyte Typing V. New York: Oxford University Press; 1995. (Biology)

Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. New York: Oxford University Press; 1995. (Biology) Tryphonas H, Lacroix F, Hayward S, Izaguirre C, Parenteau M, Fournier J. Cell surface marker evaluation of infant Macaca monkey leukocytes in peripheral whole blood using simultaneous dual-color immunophenotypic analysis. *J Med Primatol.* 1996; 25(2):89-105. (Biology)

Verdier F, Aujoulat M, Condevaux F, Descotes J. Determination of lymphocyte subsets and cytokine levels in cynomolgus monkeys. *Toxicology*. 1995; 105(1):81-90. (Biology)

Wilson AD, Shooshtari M, Finerty S, Watkins P, Morgan AJ. Selection of monoclonal antibodies for the identification of lymphocyte surface antigens in the New World primate Saguinus oedipus (cotton top tamarin). *J Immunol Methods*. 1995; 178(2):195-200. (Biology)

560811 Rev. 1 Page 2 of 2