

Technical Data Sheet

Rat IL-2 Flex Set

Product Information

Material Number:561420Size:100 TestsBead Position:E5

Assay Range: 40-10,000 pg/mL
Reactivity: QC Testing: Rat

Component Description: Rat IL-2 Standard Component Mat. No: 51-9007072

Component Storage Buffer: Lyophilized in an aqueous buffered solution containing BSA

and ProClin™ 150.

Component Description: Rat IL-2 PE Detection Reagent

Component Mat. No: 51-9007078

Component Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09%

sodium azide.

Component Description: Rat IL-2 Capture Bead E5

Component Mat. No: 51-9007080

Component Storage Buffer: Aqueous buffered solution containing fetal bovine serum

and ≤0.09% sodium azide.

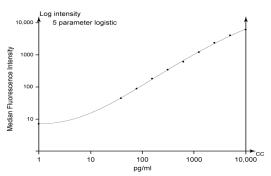


Figure 1. Example BD CBA Rat IL-2 Flex Set standard curve. Data acquired on a BD FACSArray bioanalyzer and analyzed using FCAP Array Software.

Description

The BDTM CBA Rat IL-2 Flex Set is a bead-based immunoassay capable of measuring rat interleukin-2 (IL-2) in serum and cell culture supernatant samples. Rat reactivity was determined by testing samples with the BD CBA Rat IL-2 Flex Set. The biology and function of IL-2 has been extensively reviewed in the literature. For more information on bead-based immunoassays, refer to the product insert for the BD CBA Mouse/Rat Soluble Protein Master Buffer Kit (Cat. No. 558266, 100 tests, or 558267, 500 tests).

Preparation and Storage

This BD™ CBA Flex Set contains one vial each of Capture Bead and PE Detection Reagent and two vials of Standard. The Capture Bead and PE Detection Reagent components of this flex set have been formulated to a 50x concentration to ensure product performance when multiplexed. The Standard component is lyophilized and should be transferred to a 15 mL polypropylene tube for reconstitution. When reconstituted in 4.0 mL Assay Diluent, the standard has a protein concentration of 10,000 pg/mL. Discard unused reconstituted standard, do not store or reuse. Store lyophilized standard and other components at 4°C. Protect Capture Beads and the PE Detection Reagent from prolonged exposure to light.

Application Notes

Recommended Assay Procedure: The BD CBA Rat IL-2 Flex Set must be used in conjunction with a BD CBA Mouse/Rat Soluble Protein Master Buffer Kit (Cat. No. 558266, 100 tests, or 558267, 500 tests), a flow cytometer, and FCAP Array™ Software. Detailed instructions on the use of this product can be found in the manual for the BD CBA Mouse/Rat Soluble Protein Master Buffer Kit. When following the directions in the Master Buffer Kit, the top standard point for the BD CBA Rat IL-2 Flex Set will be 10,000 pg/mL. An example standard curve is shown in figure 1.

The BD CBA Rat IL-2 Flex Set should not be used in the same assay well with any non-BD CBA Rat Soluble Protein Flex Set reagents (such as

BD Biosciences

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 800.268.5430
 32.2.400.98.95
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale. Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



561420 Rev. 3 Page 1 of 2

BD CBA Human or Mouse Soluble Protein or Cell Signaling Flex Sets). For an updated assay compatibility chart for the BD CBA Rat Soluble Protein Flex Sets, please refer to the BD CBA Flex Set System homepage at http://www.bdbiosciences.com/cbasetup.

Performance

Limit of Detection: The theoretical limit of detection is 0.46 pg/mL and was determined by evaluating the estimated result of the average MFI of the negative control (0 pg/mL, n=30) + 2 standard deviations.

Specificity		Inter-Assay Reproducibility			Intra-Assay Reproducibility		
		Mean (pg/ml)	Standard Deviation	%CV	Mean (pg/ml)	Standard Deviation	%CV
Rat IL-2	Sample 1	134.4	12.8	10%	118.4	4.9	4%
	Sample 2	571.7	30.7	5%	548.6	26.9	5%
	Sample 3	2433.1	176.1	7%	2338.5	221.9	9%

Reproducibility: The intra-assay and inter-assay reproducibility were determined for the BD CBA Rat IL-2 Flex Set by evaluating ten replicates of three different sample levels (intra-assay) and three replicates of three different sample levels from four separate experiments (inter-assay).

		ulture natant	Serum		
Sample Dilution	Average % Recovery	Range	Average % Recovery	Range	
Rat IL-2	90%	85 - 96%	71%	63 - 75%	

Recovery: Cell culture supernatant and pooled rat serum were spiked with three different levels of protein. The spiked samples were assayed and the results were compared with expected values. Serum samples were diluted 1:4 before the protein was spiked into each. Serum used was commercially available pooled rat serum.

Committee Dilection		ulture natant	Serum		
Sample Dilution	Detected (pg/ml)	% of Expected	Detected (pg/ml)	% of Expected	
Spiked sample	4521.5	100%	3461.7	100%	
1:2	2376.6	105%	2103.6	122%	
1:4	1082.2	96%	1033.2	119%	

Linearity: Cell culture supernatant and 1:4 diluted pooled rat serum were spiked with protein and serially diluted. The diluted samples were assayed and the results were compared with the original spiked sample.

Product Notices

- 1. ProClin is a trademark of Rohm and Haas Company.
- Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 4. Warning: CBA lyophilized standard contains 0.02% (w/w) of a CMIT/MIT mixture (3:1), which is a mixture of: 5-chloro-2-methyl-4-isothiazolin-3-one [EC No 247-500-7] and 2-methyl-4-isothiazolin-3-one [EC No 220-239-6] (3:1). Hazard statement: May cause an allergic skin reaction.

Precautionary statements: Wear protective gloves/eye protection. Wear protective clothing. Avoid breathing mist/vapours/spray. If skin irritation or rash occurs: Get medical advice/attention. IF ON SKIN: Wash with plenty of water. Dispose of contents/container in accordance with local/regional/national/international regulations.

BD Biosciences

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbea

 877.232.8995
 800.268.5430
 32.2.400.98.95
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale. Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



561420 Rev. 3 Page 2 of 2