

Technical Data Sheet

Alexa Fluor® 647 Mouse Anti-ERK1/2 (pT202/pY204)

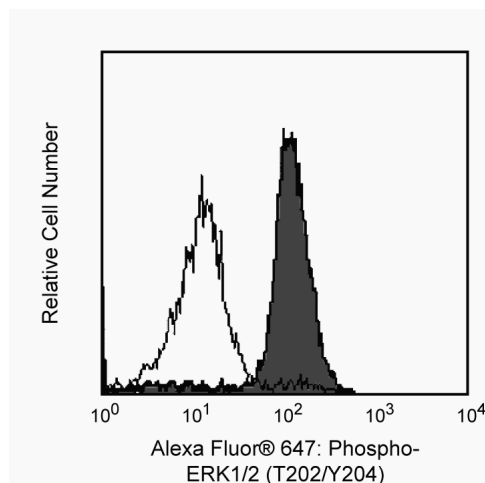
Product Information

Material Number:	561992
Alternate Name:	p44/42 MAPK; Extracellular signal-Regulated Kinase 1/2 (pT202/Y204)
Size:	250 tests
Vol. per Test:	5 µl
Clone:	20A
Immunogen:	Phosphorylated Rat ERK1 (T202/Y204) Peptide
Isotype:	Mouse IgG1
Reactivity:	QC Testing: Human Tested in Development by Western blot with the purified antibody (Cat. No. 612358): Mouse, Rat
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The members of the Mitogen-Activated Protein Kinase (MAPK) family are components of a key signal transduction cascade that links events at the cell surface to responses in the nucleus. The signaling cascade is found in species as varied as yeast and humans, with many of the proteins being well conserved. In mammals the most widely studied members of the cascade are the Extracellular signal-Regulated Kinases, ERK1 (p44 MAPK) and ERK2 (p42 MAPK). ERK1 and ERK2 share 85% homology and are activated by extracellular signals such as growth factors, hormones, and phorbol esters. Activation occurs through a series of phosphorylations by kinases activating other kinases and eventually leading to phosphorylation of the ERKs. Growth factor stimulation leads to activation of Ras and Raf, leading to phosphorylation of MEK1 (MAPK/ERK kinase) which, in turn, activates the ERKs via dual phosphorylation. Once activated, the ERKs phosphorylate other cytoplasmic signalling molecules, cell-surface receptors, microtubule-associated proteins, and transcription factors in the nucleus. Thus, the active ERK has myriad downstream effectors that implicate it in the control of cell proliferation and differentiation, as well as regulation of the cytoskeleton. Furthermore, studies have shown that elevated ERK activity is associated with some cancers.

The 20A monoclonal antibody recognizes the phosphorylated threonine 202 and tyrosine 204 (pT202/pY204) of human ERK1 and pT184/pY186 of human ERK2. The orthologous phosphorylation sites in murine ERK1 and ERK2 are T203/Y205 and T183/Y185.

**Flow cytometric analysis of ERK1/2 (pT202/pY204).**

Human peripheral blood mononuclear cells (PBMC) were either left unstimulated (unshaded) or stimulated (shaded) with 40 nM PMA for 10 minutes at 37°C. Cells were fixed with BD Cytofix™ buffer (Cat. No. 554655) for 10 minutes at 37°C and then permeabilized by adding BD Phosflow™ Perm Buffer III (Cat. No. 558050) for 30 minutes on ice. Cells were then washed twice in BD Pharmingen™ Stain Buffer (Cat. No. 554656) and stained with the Alexa Fluor® 647 mouse anti-ERK1/2 (pT202/pY204) antibody. Cells were analyzed on a BD FACSCalibur™ flow cytometry instrument. For intracellular staining of human whole blood, BD™ Phosflow Lyse/Fix buffer (Cat. No. 558049) may be used for fixation.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Application Notes

Application

Intracellular staining (flow cytometry)

Routinely Tested

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.268.5430	32.53.720.550	0120.8555.90	65.6861.0633	0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



Recommended Assay Procedure:

For more information about BD Phosflow™: Please refer to <http://www.bdbiosciences.com/support/resources/phosflow/index.jsp>

Investigators may also find the following protocols to be helpful:

Phosflow protocol for human PBMC: Please refer to http://www.bdbiosciences.com/documents/Phosflow_Protocol_for_Human_PBMCs.pdf

Phosflow protocol for human whole blood: Please refer to

http://www.bdbiosciences.com/documents/Phosflow_Protocol_for_Human_Whole_Blood_Samples.pdf

Phosflow protocol for adherent cells: Please refer to http://www.bdbiosciences.com/support/resources/protocols/protocol_adherent.jsp

Suggested Companion Products

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 ml	(none)
554655	Fixation Buffer	100 ml	(none)
557870	Fix Buffer I	250 ml	(none)
558049	Lyse/Fix Buffer 5X	250 ml	(none)
558050	Perm Buffer III	125 ml	(none)
612358	Purified Mouse Anti-ERK1/2 (pT202/pY204)	50 µg	20A
557714	Alexa Fluor® 647 Mouse IgG1 κ Isotype Control	100 tests	MOPC-21

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100-µl experimental sample (a test).
2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
3. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
5. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
7. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
8. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.

References

Boulton TG, Cobb MH. Identification of multiple extracellular signal-regulated kinases (ERKs) with antipeptide antibodies. *Cell Regul.* 1991; 2(5):357-371. (Biology)
Clark EA, Hynes RO. Ras activation is necessary for integrin-mediated activation of extracellular signal-regulated kinase 2 and cytosolic phospholipase A2 but not for cytoskeletal organization. *J Biol Chem.* 1996; 271(25):14814-14818. (Biology)
Sivaraman VS, Wang H, Nuovo GJ, Malbon CC. Hyperexpression of mitogen-activated protein kinase in human breast cancer. *J Clin Invest.* 1997; 99(7):1478-1483. (Biology)