# Technical Data Sheet

# Alexa Fluor® 647 Mouse Anti-Stat1 (pY701)

#### **Product Information**

**Material Number:** 562070 Size: 250 tests 5 µl Vol. per Test: Clone:

Immunogen: Phosphorylated Human Stat1 Peptide

**Isotype:** Mouse IgG2a Reactivity: QC Testing: Human


Tested in Development: Mouse

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

#### Description

Stat (Signal transducer and activators of transcription) proteins are critical mediators of the biologic activity of cytokines, including interleukins, interferons, erythropoietin, and growth factors. Ligand-receptor interaction leads to activation of constitutively associated JAK family kinases and subsequent recruitment/activation of Stat proteins by tyrosine phosphorylation. Active Stat proteins then move to the nucleus to promote transcription of cytokine-inducible genes. Seven Stat proteins have been cloned, each of which is differentially expressed and/or activated in a cytokine-specific and cell type-specific manner. Stat1 and Stat2 are components of the ISGF3 (Interferon-Stimulated Gene Factor 3) complex, which is the primary transcription activator induced by the binding of the interferon to a specific cell-surface receptor. Stat1 has two alternatively spliced isoforms, 91-kDa Stat1α and 84-kDa Stat1β; Stat1α has 38 additional C-terminal amino acids. In response to the binding of IFNα, IFNγ, EGF, PDGF, or CSF-1 to their respective receptors, the Stat1 subunits become tyrosine-phosphorylated at Y701, and the complex is translocated to the nucleus. This results in the formation of an active complex that includes the DNA-binding p48 subunit. This complex is responsible for modulating the transcription of the interferon-stimulated genes

The 4a monoclonal antibody recognizes the phosphorylated Y701 in Stat1α and Stat1β.

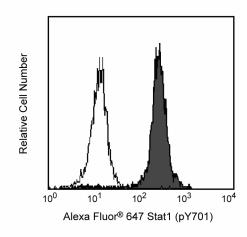

Flow cytometric analysis of Stat1 (pY701). U-937 cells (Human histiocytic lymphoma; ATCC CRL-1593.2) were either left unstimulated (unshaded) or stimulated (shaded) with 1000 U/mL recombinant human IFN-y (Cat. No. 554617) for 15 minutes at 37°C. Cells were fixed with BD Cytofix<sup>™</sup> buffer (Cat. No. 554655) for 10 minutes at 37°C and then permeabilized by adding BD Phosflow  $^{\text{\tiny TM}}$  Perm Buffer III (Cat. No. 558050 for 30 minutes on ice. Cells were then washed twice in BD Pharmingen  $^{\text{TM}}$  Stain Buffer (Cat. No. 554656) and stained with the Alexa Fluor® 647 mouse anti-Stat1 (pY701) antibody. Cells were analyzed on a BD  $FACSCalibur^{ ext{TM}}$  flow cytometery instrument. For intracellular staining of human whole blood, BD Phosflow™ Lyse/Fix buffer (Cat. No. 558049) may be used for fixation.

### **Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

#### **Application Notes**

Application

Intracellular staining (flow cytometry)

Routinely Tested

# **BD Biosciences**

bdbiosciences.com

**United States** Asia Pacific Europe Japan 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how\_to\_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



562070 Rev. 1 Page 1 of 2

#### **Recommended Assay Procedure:**

For more information about BD Phosflow<sup>TM</sup>: Please refer to http://www.bdbiosciences.com/support/resources/phosflow/index.jsp Investigators may also find the following protocols to be helpful:

Phosflow protocol for human PBMC: Please refer to http://www.bdbiosciences.com/documents/Phosflow\_Protocol\_for\_Human\_PBMCs.pdf Phosflow protocol for human whole blood: Please refer to

http://www.bdbiosciences.com/documents/Phosflow Protocol for Human Whole Blood Samples.pdf

Phosflow protocol for adherent cells: Please refer to http://www.bdbiosciences.com/support/resources/protocols/protocol\_adherent.jsp

## **Suggested Companion Products**

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 ml	(none)
554655	Fixation Buffer	100 ml	(none)
558050	Perm Buffer III	125 ml	(none)
558049	Lyse/Fix Buffer 5X	250 ml	(none)

## **Product Notices**

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10<sup>6</sup> cells in a 100-µl experimental sample (a test).
- 2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 3. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
- 4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 5. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
- 6. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 8. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

#### References

Bromberg J, Darnell JE. The role of STATs in transcriptional control and their impact on cellular function. *Oncogene*. 2000; 19(21):2468-2473. (Biology) Darnell JE Jr. STATs and gene regulation. *Science*. 1997; 277(5332):1630-1635. (Biology)

Fu XY, Zhang JJ. Transcription factor p91 interacts with the epidermal growth factor receptor and mediates activation of the c-fos gene promoter. *Cell.* 1993; 74(6):1135-1145. (Biology)

562070 Rev. 1 Page 2 of 2