

Technical Data Sheet


PE Mouse anti-Human GARP

Product Information

Material Number:	562150
Alternate Name:	LRRC32; Garpin; Glycoprotein A repetitions predominant
Size:	100 tests
Vol. per Test:	5 µl
Clone:	7B11 (also known as CMSSC-7B11)
Isotype:	Mouse IgG2b, κ
Reactivity:	QC Testing: Human
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The 7B11 (also known as CMSSC-7B11) monoclonal antibody specifically binds to human GARP (Glycoprotein A repetitions predominant). The LRRC32 (Leucine rich repeat containing 32) gene encodes the 662 amino acid-residue, 80 kDa transmembrane GARP glycoprotein that has an extracellular region composed primarily of 20 leucine-rich repeats. GARP is specifically expressed on Treg cells activated through the T cell receptor (TCR). Ectopic expression of GARP in human naïve T cells inhibited their proliferation and cytokine secretion upon TCR activation. Remarkably, GARP over-expression in naïve T cells induced expression of FoxP3 and endowed them with a partial suppressive function. The extracellular, but not the cytoplasmic region, of GARP was necessary for these functions. GARP serves as a receptor for latent TGF-beta which may play a role in the suppressive action of Treg cells. GARP is also expressed on platelets and other tissues, however the function on these cells is not known.

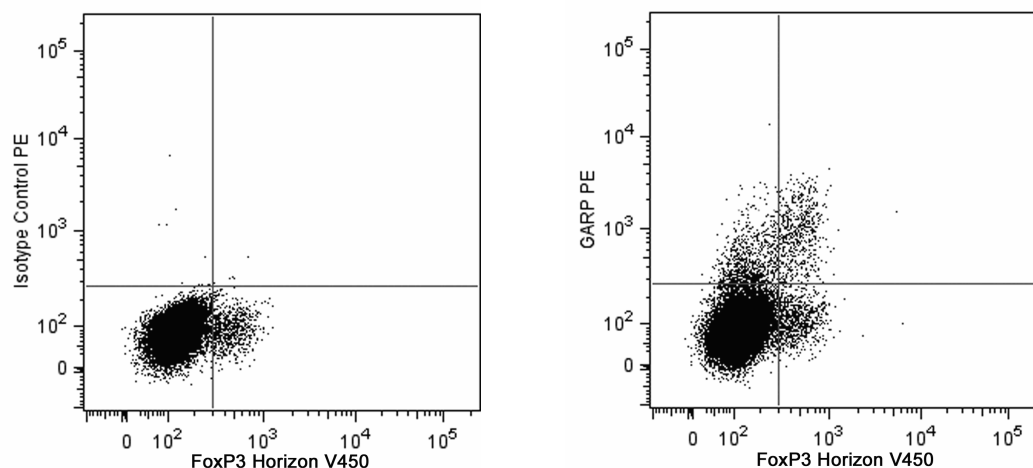

Multicolor flow cytometric analysis of GARP expressed on activated human peripheral blood lymphocytes. Human peripheral blood mononuclear cells (PBMC) were activated (12 hours) with plate-bound Purified NA/LE Mouse Anti-Human CD3 (Cat. No. 555329) and Purified NA/LE Mouse Anti-Human CD28 (Cat. No. 555725). The activated PBMC were washed and stained with either PE Mouse IgG2b, κ Isotype Control (Cat. No. 556656) or PE Mouse Anti-Human GARP (Cat. No. 562150). The cells were then washed and stained with BD Horizon™ V450 Mouse Anti-Human FoxP3 (Cat. No. 560459) using FoxP3 staining protocol. Two-color flow cytometric dot plots showing the correlated expression staining patterns of V450 FoxP3 versus PE anti-Human GARP (right panel) or PE Mouse IgG2b, κ isotype control (left panel) were derived from events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed.

Application Notes

Application

Flow cytometry

Routinely Tested

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.268.5430	32.53.720.550	0120.8555.90	65.6861.0633	0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



Recommended Assay Procedure:**Suggested Staining Procedure for PE Mouse anti-Human GARP Antibody:**

1. Harvest PBMC after stimulation (12 hours) with plate-bound Purified NA/LE Mouse Anti-Human CD3 (Cat. No. 555329) and Purified NA/LE Mouse Anti-Human CD28 (Cat. No. 555725).
 2. Wash the cells twice with stain buffer (eg. BD Pharmingen™ Stain Buffer (FBS), Cat. No. 554656).
 3. Stain 1×10^6 cells either with the PE Mouse anti-Human GARP antibody (Cat. No. 562150) or with PE Mouse IgG2b, κ Isotype control (Cat. No. 556656) for 30 minutes on ice, protected from light.
 4. Wash cells twice with stain buffer.
 5. Stain with BD Horizon™ V450 Mouse anti-Human FoxP3. Refer to the Technical Data Sheet of Cat. No. 560459 for detailed FoxP3 staining protocol.
- In brief,
- a. Add 2 ml of $1 \times$ FoxP3 buffer A to the cell pellet.
 - b. Centrifuge and incubate in 0.5 ml of buffer C for 30 minutes.
 - c. Wash twice with stain buffer and stain with V450 FoxP3 antibody (Cat. No. 560459) for 30-45 minutes.
 - d. Wash twice with stain buffer and acquire on the Flow cytometer.

Suggested Companion Products

Catalog Number	Name	Size	Clone
556656	PE Mouse IgG2b, κ Isotype Control	50 tests	27-35
554656	Stain Buffer (FBS)	500 ml	(none)
555329	Purified NA/LE Mouse Anti-Human CD3	0.5 mg	UCHT1
555725	Purified NA/LE Mouse Anti-Human CD28	0.5 mg	CD28.2
560459	V450 Mouse Anti-Human FoxP3	120 tests	259D/C7
552843	Anti-Mouse Ig, κ /Negative Control (FBS) Compensation Particles Set	6.0 ml	187.1
555899	Lysing Buffer	100 ml	(none)
560098	Human FoxP3 Buffer Set	100 tests	(none)

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
4. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
6. BD Horizon™ V450 has a maximum absorption of 406 nm and maximum emission of 450 nm. Before staining with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
7. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

Ollendorff V, Noguchi T, deLapeyriere O, Birnbaum D. The GARP gene encodes a new member of the family of leucine-rich repeat-containing proteins. *Cell Growth Differ.* 1994; 5(2):213-219. (Biology)

Stockis J, Colau D, Coulie PG, Lucas S. Membrane protein GARP is a receptor for latent TGF-beta on the surface of activated human Treg. *Eur J Immunol.* 2009; 39(12):3315-3322. (Biology)