Technical Data Sheet

PE-CF594 Mouse Anti-Human CD20

Product Information

Material Number: 562550

Alternate Name: MS4A1; B1; Bp35; LEU-16; S7

 Size:
 50 tests

 Vol. per Test:
 5 μl

 Clone:
 2H7

Isotype: Mouse (BALB/c) IgG2b, κ

Reactivity: Human

QC Testing: Rhesus or Cynomolgus Macaque or Baboon


Workshop: IV B201

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The 2H7 monoclonal antibody specifically binds to CD20 that is encoded by the MS4A1 (Membrane-spanning 4-domains, subfamily A, member 1) gene. CD20 is a 33-37 kDa unglycosylated four-transmembrane phosphoprotein. CD20 is expressed on pre-B-cells, resting and activated B cells and follicular dendritic cells but not on plasma cells. Low level CD20 expression is observed on a small subset of normal circulating T lymphocytes. The CD20 molecule is involved in the regulation of B-cell activation.

This antibody is conjugated to BD HorizonTM PE-CF594, which has been developed exclusively by BD Biosciences as a better alternative to PE-Texas Red®. PE-CF594 excites and emits at similar wavelengths to PE-Texas Red® yet exhibits improved brightness and spectral characteristics. Due to PE having maximal absorption peaks at 496 nm and 564 nm, PE-CF594 can be excited by the blue (488-nm), green (532-nm) and yellow-green (561-nm) lasers and can be detected with the same filter set as PE-Texas Red® (eg 610/20-nm filter).

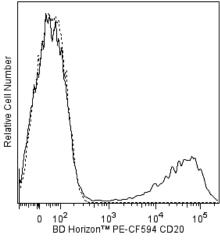

Flow cytometric analysis of CD20 expression on Rhesus macaque peripheral blood lymphocytes. Rhesus macaque whole blood was stained with BD Horizon™ PE-CF594 Mouse Anti-Human CD20 antibody (Cat. No. 562550; solid line histogram) or with a BD Horizon™ PE-CF594 Mouse IgG2b, κ Isotype Control (Cat. No. 562305; dashed line histogram). The erythrocytes were lysed with BD Pharm Lyse™ Lysing Buffer (Cat. No. 555899). The fluorescence histograms were derived from events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The antibody was conjugated with BD Horizon™ PE-CF594 under optimum conditions, and unconjugated antibody and free PE-CF594 were removed.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number	Name	Size	Clone	
562305	PE-CF594 Mouse IgG2b, κ Isotype Control	0.1 mg	27-35	
555899	Lysing Buffer	100 ml	(none)	

BD Biosciences

bdbiosciences.com

United States Canada Europe Japan Asia Pacific Latin America/Caribbean 877.232.8995 800.979.9408 32.53.720.550 0120.8555.90 65.6861.0633 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



Product Notices

- Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use $1 \times 10^{\circ}6$ cells in a 100- μ l experimental
- 3. An isotype control should be used at the same concentration as the antibody of interest.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- Please observe the following precautions: Absorption of visible light can significantly alter the energy transfer occurring in any tandem fluorochrome conjugate; therefore, we recommend that special precautions be taken (such as wrapping vials, tubes, or racks in aluminum foil) to prevent exposure of conjugated reagents, including cells stained with those reagents, to room illumination.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
- Texas Red is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- CFTM is a trademark of Biotium, Inc.
- When excited by the yellow-green (561-nm) laser, the fluorescence may be brighter than when excited by the blue (488-nm) laser.
- 11. This product is provided under an Agreement between BIOTIUM and BD Biosciences. The manufacture, use, sale, offer for sale, or import of this product is subject to one or more patents or pending applications owned or licensed by Biotium, Inc. This product, and only in the amount purchased by buyer, may be used solely for buyer's own internal research, in a manner consistent with the accompanying product literature. No other right to use, sell or otherwise transfer (a) this product, or (b) its components is hereby granted expressly, by implication or by estoppel. This product is for research use only. Diagnostic uses require a separate license from Biotium, Inc. For information on purchasing a license to this product including for purposes other than research, contact Biotium, Inc., 3159 Corporate Place, Hayward, CA 94545, Tel: (510) 265-1027. Fax: (510) 265-1352. Email: btinfo@biotium.com.
- 12. Because of the broad absorption spectrum of the tandem fluorochrome, extra care must be taken when using multi-laser cytometers, which may directly excite both PE and CFTM594.

References

Hultin LE, Hausner MA, Hultin PM, Giorgi JV. CD20 (pan-B cell) antigen is expressed at a low level on a subpopulation of human T lymphocytes. Cytometry. 1993; 14(2):193-204. (Biology)

Knapp W, Dorken B, Rieber EP, et al, ed. Leucocyte Typing IV. New York: Oxford University Press; 1989:1-1208. (Clone-specific) Schlossman S, Boumell L, et al, ed. Leucocyte Typing V. New York: Oxford University Press; 1995. (Biology)

BD Biosciences

bdbiosciences.com

United States Canada Europe Asia Pacific Latin America/Caribbean 800.979.9408 32.53.720.550 0120.8555.90 65.6861.0633 877.232.8995 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation

of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited.
For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



562550 Rev. 1 Page 2 of 2