

Technical Data Sheet


BV421 Mouse Anti-Human CD34**Product Information**

Material Number:	562577
Alternate Name:	gp105-120; My10; Hematopoietic progenitor cell antigen CD34
Size:	100 tests
Vol. per Test:	5 µl
Clone:	581
Isotype:	Mouse IgG1, κ
Reactivity:	QC Testing: Human
Workshop:	V MA27, VI E004
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The 581 monoclonal antibody specifically binds to CD34, a sialomucin-like type I transmembrane glycoprotein. This single-chain, 105-120 kDa, heavily O-glycosylated protein is expressed on hematopoietic progenitor cells, vascular endothelium, bone marrow stromal cells and embryonic fibroblasts. The cytoplasmic region of the CD34 antigen is a target for phosphorylation by activated protein kinase C suggesting CD34 may play a role in signal transduction. CD34 may also play a role as an adhesion molecule since it binds to CD62E and CD62L. Clone 581 binds to the class III CD34 epitope, it is resistant to neuraminidase, chymopapain and glycoprotease. The 581 antibody blocks reactivity of another anti-CD34 monoclonal antibody, 8G12.

The antibody was conjugated to BD Horizon™ BV421 which is part of the BD Horizon™ Brilliant Violet™ family of dyes. With an Ex Max of 407-nm and Em Max at 421-nm, BD Horizon™ BV421 can be excited by the violet laser and detected in the standard Pacific Blue™ filter set (eg, 450/50-nm filter). BD Horizon™ BV421 conjugates are very bright, often exhibiting a 10 fold improvement in brightness compared to Pacific Blue™ conjugates.

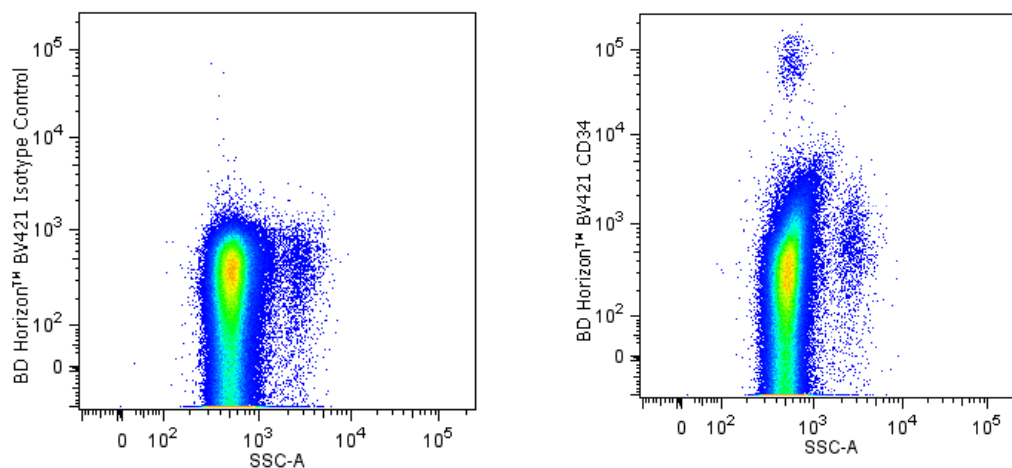

Multicolor flow cytometric analysis of CD34 expressed by human peripheral blood mononuclear cells. Human peripheral blood mononuclear cells were stained with APC Mouse Anti-Human CD14 (Cat. No. 555399/561708) and either a BD Horizon™ BV421 Mouse IgG1, Isotype Control (Cat No. 562438, Left Panel) or with the BD Horizon™ BV421 Mouse Anti-Human CD34 antibody (Cat. No. 562577, Right Panel). Flow cytometric dot plots showing side-scattered light versus CD34 (or Ig isotype control staining) were derived from gated events based on light scattering characteristics for CD14-negative cells. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ BV421 under optimum conditions, and unconjugated antibody and free BD Horizon™ BV421 were removed.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



Application Notes

Application

Flow cytometry

Routinely Tested

Suggested Companion Products

Catalog Number	Name	Size	Clone
562438	BV421 Mouse IgG1, k Isotype Control	50 µg	X40
554656	Stain Buffer (FBS)	500 ml	(none)
555399	APC Mouse Anti-Human CD14	100 tests	M5E2
561708	APC Mouse Anti-Human CD14	25 tests	M5E2

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100-µl experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
4. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.
5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
7. Pacific Blue™ is a trademark of Molecular Probes, Inc., Eugene, OR.
8. Brilliant Violet™ 421 is a trademark of Sirigen.

References

Egeland T, Tjonnfjord G, Steen R, Gaudernack G, Thorsby E. Positive selection of bone marrow-derived CD34 positive cells for possible stem cell transplantation. *Transplant Proc.* 1993; 25(1):1261-1263. (Biology)

Knapp W, Dörken B, Gilks WR, et al, ed. *Leukocyte Typing IV*. New York, NY: Oxford University Press; 1989:1-1182. (Biology)

Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. New York: Oxford University Press; 1995. (Clone-specific)

Zola H, Swart B, Nicholson I, Voss E. *Leukocyte and Stromal Cell Molecules. The CD Markers*. Hoboken, New Jersey: John Wiley & Sons, Inc.; 2007:1-581. (Biology)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

