Technical Data Sheet

Alexa Fluor® 647 Mouse Anti-Human CD107A

Product Information

Material Number: 562622

Alternate Name: LAMP1; LAMP-1; LAMPA; LGP120


Immunogen: Human Adult Adherent Peripheral Blood Cells

Workshop: V P008

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The H4A3 monoclonal antibody specifically binds to the heavily glycosylated 110 kDa Lysosomal-associated membrane protein 1, LAMP-1. LAMP-1 is a widely expressed intracellular antigen. It is also expressed on the surface of activated platelets, PHA-activated lymphocytes, cytotoxic T cells and NK cells, and some tumor cell lines, including U937 and KG1a. LAMP-1 has been shown to be a ligand for E-selectin-mediated cell adhesion. LAMP-1 and LAMP-2 (CD107b) are carriers for poly-N-acetyllactosamines and are able to display sialyl Le[x] termini.

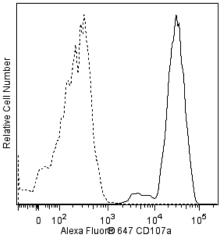

Flow cytometric analysis of CD107a expressed by Jurkat cells. Jurkat cells were fixed with BD Cytofix™ Fixation Buffer (Cat. No. 554655) and permeabilized with BD Perm/Wash™ Buffer (Cat. No. 554723) and subsequently stained either with Alexa Fluor® 647 Mouse [gG1, κ Isotype Control (Cat. No. 557732; dashed line histogram) or with the Alexa Fluor® 647 Mouse Anti-Human CD107a antibody (Cat. No. 562622; solid line histogram). The fluorescence histograms were derived from events with the forward and side light-scatter characteristics of intact Jurkat cells. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Application Notes

Application

Intracellular staining (flow cytometry)

Routinely Tested

Suggested Companion Products

Catalog Number	Name Name	Size	Clone
557732	Alexa Fluor® 647 Mouse IgG1 κ Isotype Control	100 tests	MOPC-21
554656	Stain Buffer (FBS)	500 ml	(none)
554655	Fixation Buffer	100 ml	(none)
554723	Perm/Wash Buffer	100 ml	(none)

BD Biosciences

bdbiosciences.com

United States Canada Europe Japan Asia Pacific Latin America/Caribbean 877.232.8995 800.979.9408 32.53.720.550 0120.8555.90 65.6861.0633 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



562622 Rev. 1 Page 1 of 2

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁻⁶ cells in a 100-μl experimental sample (a test).
- An isotype control should be used at the same concentration as the antibody of interest. 2
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols. 3.
- The Alexa Fluor®, Pacific BlueTM, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific BlueTM dye, and Cascade Blue® dye are covered by pending and issued patents.
- 5. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
- Source of all serum proteins is from USDA inspected abattoirs located in the United States.

References

Alto NM, Soderling J, Scott JD. Rab32 is an A-kinase anchoring protein and participates in mitochondrial dynamics. J Biol Chem. 2002; 158(4):659-668. (Biology) Febbraio M, Silverstein RL. Identification and characterization of LAMP-1 as an activation-dependent platelet surface glycoprotein. J Biol Chem. 1990; 265(30):18531-18537. (Biology)

Fukuda M. Lysosomal membrane glycoproteins. Structure, biosynthesis, and intracellular trafficking. J Biol Chem. 1991; 266(32):21327-21330. (Biology) Fukuda M, Viitala J, Matteson J, Carlsson SR. Cloning of cDNAs encoding human lysosomal membrane glycoproteins, h-lamp-1 and h-lamp-2. Comparison of their deduced amino acid sequences. J Biol Chem. 1988; 263(35):18920-18928. (Biology)

Hocking DC, Kowalski K. A cryptic fragment from fibronectin's III1 module localizes to lipid rafts and stimulates cell growth and contractility. J Cell Biol. 2002; 158(1):175-184. (Biology)

Mane SM, Marzella L, Bainton DF, Holt VK, Cha Y, Hildreth JE, August JT. Purification and characterization of human lysosomal membrane glycoproteins. Arch Biochem Biophys. 1989'; 268(1):360-378. (Immunogen: Electron microscopy, Flow cytometry, Immunoaffinity chromatography, Immunohistochemistry Immunoprecipitation)

Sawada R, Lowe JB, Fukuda M. E-selectin-dependent adhesion efficiency of colonic carcinoma cells is increased by genetic manipulation of their cell surface lysosomal membrane glycoprotein-1 expression levels. J Biol Chem. 1993; 268(17):12675-12681. (Biology)

Schlossman SF, Boumsell L, Gilks W, et al, ed. Leukocyte Typing V: White Cell Differentiation Antigens. New York: Oxford University Press; 1995

Spoerl Z, Stumpf M, Noegel AA, Hasse A. Oligomerization, F-actin interaction, and membrane association of the ubiquitous mammalian coronin 3 are mediated by its carboxyl terminus. J Biol Chem. 2002; 277(50):48858-48867. (Biology)

BD Biosciences

bdbiosciences.com

United States Canada Europe Asia Pacific Latin America/Caribbean 800.979.9408 32.53.720.550 0120.8555.90 65.6861.0633 877.232.8995 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation

of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express

written authorization of Becton, Dickinson and Company is stictly prohibited.
For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



562622 Rev. 1 Page 2 of 2