Technical Data Sheet

FITC Mouse anti-Human CD56

Product Information

Material Number: 562794

Alternate Name: NCAM1; NCAM-1; Neural cell adhesion molecule 1; NCAM; NKH1; MSK39

 Size:
 100 tests

 Vol. per Test:
 5 μl


 Clone:
 B159

Workshop: V NK75

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The B159 monoclonal antibody specifically binds to CD56. CD56 is a heavily glycosylated adhesion protein that is present on a subpopulation of peripheral blood large granular lymphocytes that demonstrate natural killer activity. CD56 is also expressed on a subset of T cells but is not expressed on myeloid cells, erythrocytes or B cells. This antigen is a pan-NK-cell marker. CD56 is virtually identical to an isoform of the neutral cell adhesion molecule (NCAM), a structure mediating homotypic and heterotypic cell-cell interactions.

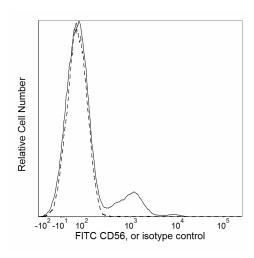

Flow cytometric analysis of CD56 expression on human peripheral blood lymphocytes. Human whole blood was stained with the FITC Mouse anti-Human CD56 antibody (Cat. No. 562794; solid line histogram) or with FITC Mouse Ig61, к Isotype Control (Cat. No. 555748; dashed line histogram). The erythrocytes were lysed with BD Pharm Lyse™ Lysing Buffer (Cat. No. 555899). The fluorescence histograms were derived from gated events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometry was performed using a BD FACSCanto™ II Flow Cytometer System.

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed. Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

| Catalog Number | Name | Size | Clone | |
|----------------|------------------------------------|-----------|---------|--|
| 555748 | FITC Mouse IgG1, κ Isotype Control | 100 tests | MOPC-21 | |
| 554656 | Stain Buffer (FBS) | 500 ml | (none) | |
| 555899 | Lysing Buffer | 100 ml | (none) | |

BD Biosciences

bdbiosciences.com

United States Canada Europe Japan Asia Pacific Latin America/Caribbean 877.232.8895 800.979.9408 32.53.720.550 0120.8555.90 65.6861.0633 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



562794 Rev. 1 Page 1 of 2

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use $1 \times 10^{\circ}6$ cells in a 100- μ l experimental sample (a test).
- An isotype control should be used at the same concentration as the antibody of interest. 2
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols. 3.
- For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
- 5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- Source of all serum proteins is from USDA inspected abattoirs located in the United States. 6.

References

Bennett IM, Zatsepina O, Zamai L, Azzoni L, Mikheeva T, Perussia B. Definition of a natural killer NKR-P1A+/CD56-/CD16- functionally immature human NK cell subset that differentiates in vitro in the presence of interleukin 12. J Exp Med. 1996; 184(5):1845-1856. (Biology)

Campbell JJ, Qin S, Unutmaz D, Soler D, Murphy KE, Hodge MR, Wu L, Butcher EC. Unique subpopulations of CD56+ NK and NK-T peripheral blood lymphocytes identified by chemokine receptor expression repertoire. J Immunol. 2001; 166(11):6477-6482. (Biology)

Cooper MA, Fehniger TA, Caligiuri MA. The biology of human natural killer-cell subsets. Trends Immunol. 2001; 22(11):633-640. (Biology)

Cunningham BA, Hemperly JJ, Murray BA, Prediger EA, Brackenbury R, Edelman GM. Neural cell adhesion molecule: structure, immunoglobulin-like domains, cell surface modulation, and alternative RNA splicing. Science. 1987; 236(4803):799-806. (Biology)

Edelman GM. Cell adhesion molecules in the regulation of animal form and tissue pattern. Annu Rev Cell Biol. 1986; 2:81-116. (Biology)

Galandrini R, Tassi I, Mattia G, Lenti L, Piccoli M, Frati L, Santoni A. SH2-containing inositol phosphatase (SHIP-1) transiently translocates to raft domains and modulates CD16-mediated cytotoxicity in human NK cells. Blood. 2001; 100(13):4581-4589. (Biology)

Gerosa F, Baldani-Guerra B, Nisii C, Marchesini V, Carra G, Trinchieri G. Reciprocal activating interaction between natural killer cells and dendritic cells. J Exp Med. 2002; 195(3):327-333. (Biology)

Lanier LL, Chang C, Azuma M, Ruitenberg JJ, Hemperly JJ, Phillips JH. Molecular and functional analysis of human natural killer cell-associated neural cell adhesion molecule (N-CAM/CD56). J Immunol. 1991; 146(12):4421-4426. (Biology)

Lanier LL, Le AM, Civin CI, Loken MR, Phillips JH. The relationship of CD16 (Leu-11) and Leu-19 (NKH-1) antigen expression on human peripheral blood NK cells and cytotoxic T lymphocytes. J Immunol. 1986; 136(12):4480-4486. (Biology)

Lanier LL, Testi R, Bindl J, Phillips JH. Identity of Leu-19 (CD56) leukocyte differentiation antigen and neural cell adhesion molecule. J Exp Med. 1989; 169(6):2233-2238. (Biology)

Nitta T, Yagita H, Sato K, Okumura K. Involvement of CD56 (NKH-1/Leu-19 antigen) as an adhesion molecule in natural killer-target cell interaction. J Exp Med. 1989; 170(5):1757-1761. (Biology)

Phillips JH, Lanier LL. Dissection of the lymphokine-activated killer phenomenon: relative contribution of peripheral blood natural killer cells and T lymphocytes to cytolysis. J Exp Med. 1986; 164(3):814-825. (Biology)

Ritz J. Trinchieri G. Lanier LL. NK-cell Antigens; Section Report. In: Schlossman SF, Bournsell L, Gilks W, ed. Leucocyte Typing V, Oxford; Oxford University Press: 1995:1367-1372. (Clone-specific)

Schubert W, Zimmermann K, Cramer M, Starzinski-Powitz A. Lymphocyte antigen Leu-19 as a molecular marker of regeneration in human skeletal muscle. Proc Natl Acad Sci U S A. 1989; 86(1):307-311. (Biology)

BD Biosciences

bdbiosciences.com

United States Canada Europe Asia Pacific Latin America/Caribbean 800.979.9408 32.53.720.550 0120.8555.90 65.6861.0633 55.11.5185.9995 877.232.8995

For country contact information, visit bdbiosciences.com/contact Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation

of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



562794 Rev. 1 Page 2 of 2