


Technical Data Sheet

FITC Mouse Anti-iNOS/NOS Type II**Product Information**

Material Number:	610330
Size:	50 µg
Concentration:	250 µg/ml
Clone:	6/iNOS/NOS Type II
Immunogen:	Mouse iNOS aa. 961-1144
Isotype:	Mouse IgG2a
Reactivity:	QC Testing: Mouse
Storage Buffer:	Aqueous buffered solution containing BSA, glycerol, and ≤0.09% sodium azide.

Description

Nitric oxide synthase (NOS), a cell-type specific enzyme, catalyzes the synthesis of nitric oxide (NO). NO is a short-lived radical that transmits cellular signals involved in vasorelaxation, neurotransmission, and cytotoxicity. In macrophages and other cell types, NOS (iNOS or macNOS) activity increases following exposure to cytokines (IFN- γ , TNF- α , and IL-1) and microbial products (lipopolysaccharide (LPS)), iNOS is activated independently of Ca²⁺/calmodulin and its level of expression is tightly controlled by several transcription factors, including NF κ B. Data indicates that TGF- β affects translation of iNOS mRNA and decreases iNOS protein stability. Normally undetectable in brain tissue, iNOS mRNA has been observed in CNS tissues of animals under experimental pathologic conditions. iNOS and nNOS share 51% amino acid homology with the greatest degree of divergence in the calmodulin binding domain.

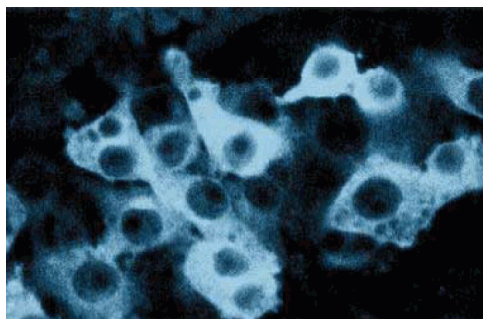

Immunofluorescence staining of iNOS/NOS Type II on mouse macrophages (RAW 264.7) treated with 10 ng/mL IFN γ + 1 µg/mL LPS for 12 hrs.

Preparation and Storage

Store undiluted at -20°C.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

Application Notes**Application**

Fluorescence microscopy

Routinely Tested

Product Notices

1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.

References

Angulo I, de las Heras FG, Garcia-Bustos JF, Gargallo D, Munoz-Fernandez MA, Fresno M. Nitric oxide-producing CD11b(+)Ly-6G(Gr-1)(+)CD31(ER-MP12)(+) cells in the spleen of cyclophosphamide-treated mice: implications for T-cell responses in immunosuppressed mice. *Blood*. 2000; 95(1):212-220. (Biology: Flow cytometry)

Koprowski H, Zheng YM, Heber-Katz E, et al. In vivo expression of inducible nitric oxide synthase in experimentally induced neurologic diseases. *Proc Natl Acad Sci U S A*. 1993; 90(7):3024-3027. (Biology)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



Moulian N, Truffault F, Gaudry-Talarmain YM, Serraf A, Berrih-Aknin S. In vivo and in vitro apoptosis of human thymocytes are associated with nitrotyrosine formation. *Blood*. 2001; 97(11):3521-3530. (Biology: Immunofluorescence)

Siedlar M, Mytar B, Krzeszowiak A, et al. Demonstration of iNOS-mRNA and iNOS in human monocytes stimulated with cancer cells in vitro. *J Leukoc Biol*. 1999; 65(5):597-604. (Biology: Flow cytometry, Immunofluorescence)

Xie QW, Cho HJ, Calaycay J, et al. Cloning and characterization of inducible nitric oxide synthase from mouse macrophages. *Science*. 1992; 256(5054):225-228. (Biology)