

CD4 antibodies, human

For research use only

One test corresponds to labeling of up to 10^7 cells in a total volume of 100 μ L.

Product	Content	Order no.
CD4-FITC	for 30 tests	130-098-160
CD4-FITC	for 100 tests	130-092-358
CD4-VioBright FITC	for 30 tests	130-104-564
CD4-VioBright FITC	for 100 tests	130-104-515
CD4-PE	for 30 tests	130-098-167
CD4-PE	for 100 tests	130-092-373
CD4-APC	for 30 tests	130-098-158
CD4-APC	for 100 tests	130-092-374
CD4-VioBlue	for 30 tests	130-098-163
CD4-VioBlue	for 100 tests	130-094-153
CD4-VioGreen	for 30 tests	130-098-168
CD4-VioGreen	for 100 tests	130-096-900
CD4-PerCP	for 30 tests	130-098-155
CD4-PerCP	for 100 tests	130-094-963
CD4-PE-Vio615	for 30 tests	130-108-305
CD4-PE-Vio615	for 100 tests	130-108-281
CD4-PE-Vio770	for 30 tests	130-098-149
CD4-PE-Vio770	for 100 tests	130-096-552
CD4-APC-Vio770	for 30 tests	130-098-153
CD4-APC-Vio770	for 100 tests	130-096-652
CD4-PerCP-Vio700	for 100 tests	130-097-580
CD4-Biotin	for 100 tests	130-098-546

Warnings

Reagents contain sodium azide. Under acidic conditions sodium azide yields hydrazoic acid, which is extremely toxic. Azide compounds should be diluted with running water before discarding. These precautions are recommended to avoid deposits in plumbing where explosive conditions may develop.

Technical data and background information

Antigen	CD4
Clone	VIT4
Isotype	mouse IgG2ak
Isotype control	Mouse IgG2a – isotype control antibodies
Alternative names of antigen	CD4mut, L3T4, Leu-3, T4
Molecular mass of antigen [kDa]	48

Distribution of antigen	dendritic cells, granulocytes, Langerhans cells, lymphocytes, macrophages, monocytes, neutrophils, T cells, T helper cells, thymocytes
Product format	Antibodies are supplied in buffer containing stabilizer and 0.05% sodium azide.
Fixation	With the exception of the VioBlue conjugate the antibody is suited for staining of formaldehyde-fixed cell. Cells should be stained prior to fixation.
Storage	Store protected from light at 2–8 °C. Do not freeze.

The monoclonal antibody VIT4 recognizes the human CD4 antigen which is highly expressed on human T helper cells and thymocytes, and at lower levels on monocytes and dendritic cells. It is responsible for the recognition of the MHC class II antigen.

The VIT4 antibodies can be used for evaluating the transfection efficiency of cells transfected with pMACS 4.1 or pMACS 4-IRES when using the MACSelect™ – 4 Transfected Cell Selection Kit.

Reagent requirements


- Buffer: Prepare a solution containing phosphate-buffered saline (PBS), pH 7.2, 0.5% bovine serum albumin (BSA), and 2 mM EDTA by diluting MACS® BSA Stock Solution (# 130-091-376) 1:20 with autoMACS® Rinsing Solution (# 130-091-222). Keep buffer cold (2–8 °C).
Note: EDTA can be replaced by other supplements such as anticoagulant citrate dextrose formula-A (ACD-A) or citrate phosphate dextrose (CPD). Buffers or media containing Ca^{2+} or Mg^{2+} are not recommended for use.
- (Optional) FcR Blocking Reagent, human (# 130-059-901) to avoid Fc receptor-mediated antibody labeling.
- (Optional) Fluorochrome-conjugated anti-biotin antibodies, e.g., Anti-Biotin-PE (# 130-090-756) as secondary antibody reagent in combination with biotinylated antibodies.
- (Optional) Propidium Iodide Solution (# 130-093-233) for flow cytometric exclusion of dead cells without fixation.
- (Optional) Fixation and Dead Cell Discrimination Kit (# 130-091-163) for cell fixation and flow cytometric exclusion of dead cells.


Protocol for cell surface staining

- The recommended antibody dilution for labeling of cells and subsequent analysis by flow cytometry is 1:11 for up to 10^7 cells/100 μL of buffer.
 - Volumes given below are for up to 10^7 nucleated cells. When working with fewer than 10^7 cells, use the same volumes as indicated. When working with higher cell numbers, scale up all reagent volumes and total volumes accordingly (e.g. for 2×10^7 nucleated cells, use twice the volume of all indicated reagent volumes and total volumes).
1. Determine cell number.
 2. Centrifuge cell suspension at 300×g for 10 minutes. Aspirate supernatant completely.
 3. Resuspend up to 10^7 nucleated cells per 100 μL of buffer.
 4. Add 10 μL of the antibody.
 5. Mix well and incubate for 10 minutes in the dark in the refrigerator (2–8 °C).
Note: Higher temperatures and/or longer incubation times may lead to non-specific cell labeling. Working on ice requires increased incubation times.
 6. Wash cells by adding 1–2 mL of buffer and centrifuge at 300×g for 10 minutes. Aspirate supernatant completely.
 7. (Optional) If biotinylated antibody was used, resuspend the cell pellet in 100 μL of buffer, add 10 μL of fluorochrome-conjugated anti-biotin antibody, and continue as described in steps 5 and 6.
 8. Resuspend cell pellet in a suitable amount of buffer for analysis by flow cytometry or fluorescence microscopy.


Examples of immunofluorescent staining

Human peripheral blood mononuclear cells (PBMCs) were stained with CD4 antibodies or with the corresponding isotype control antibodies (left image) and analyzed by flow cytometry using the MACSQuant[®] Analyzer. The Tandem Signal Enhancer has been used to increase binding specificity of tandem-dye–conjugated antibodies. Cell debris and dead cells were excluded from the analysis based on scatter signals and propidium iodide fluorescence or 4',6-diamidino-2-phenylindole (DAPI) fluorescence, as in the case of tandem conjugates.

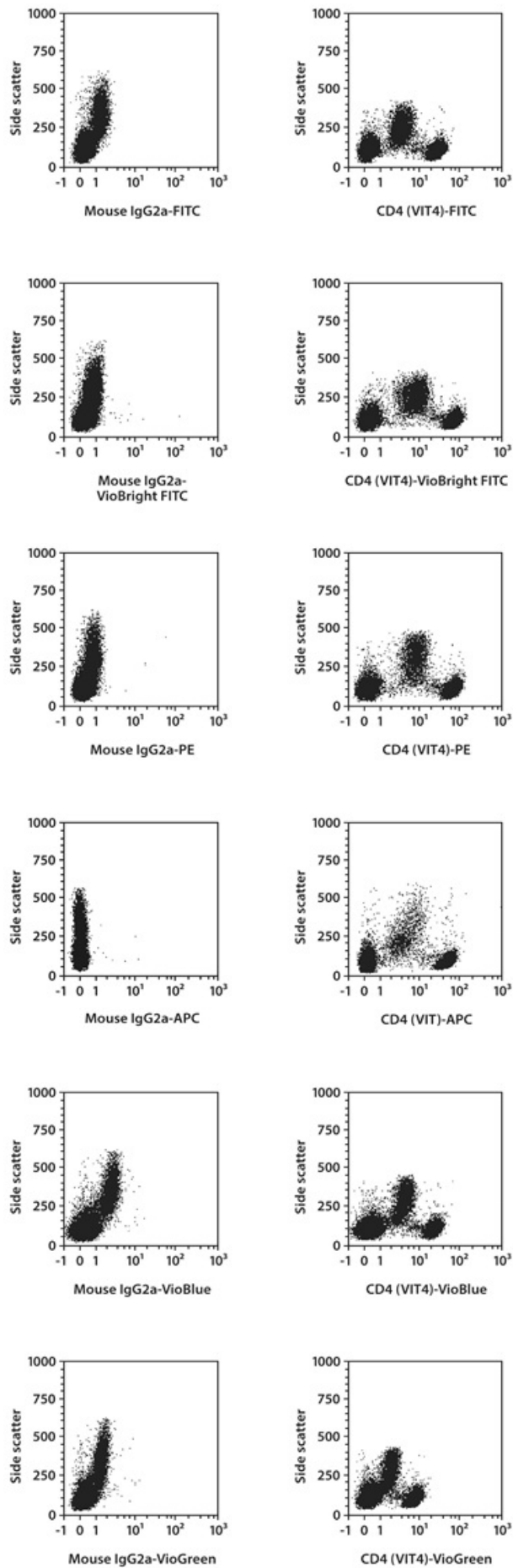

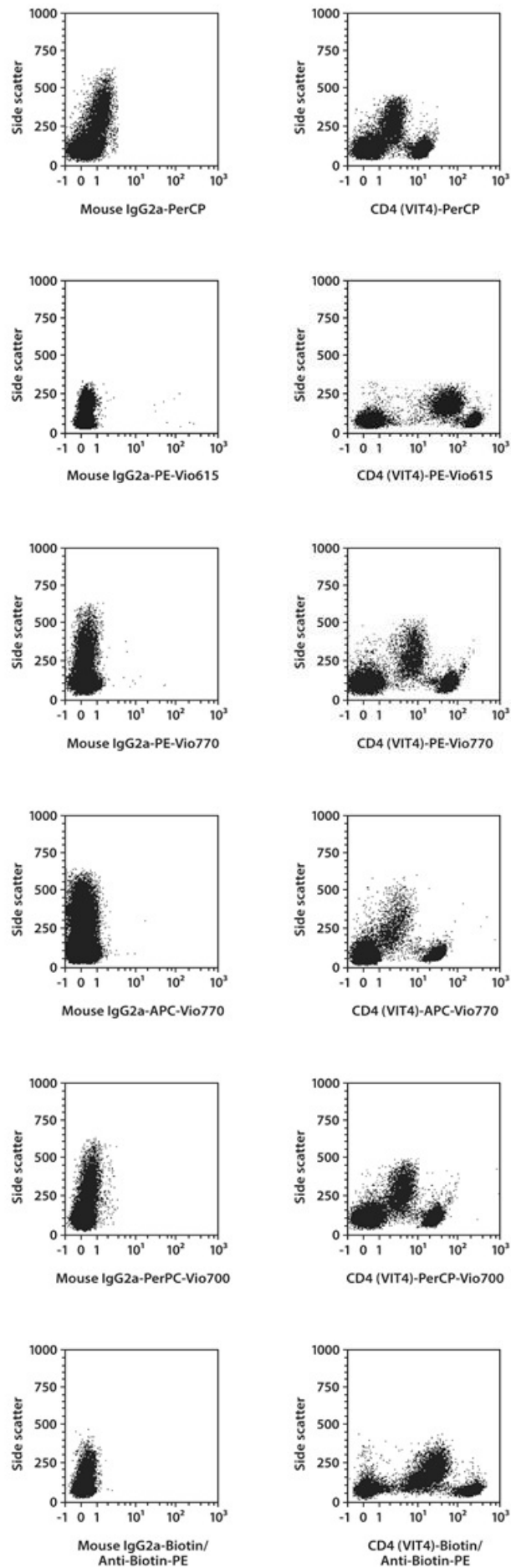

The 8 µm cryosection of the human tonsil was fixed with 100% ice-cold acetone for three minutes. The cryosection was stained with the CD4 antibodies, clone VIT4 (blue), CD279 antibodies, clone PD1.3.1.3 (red) and CD79a antibodies, clone HM47 (green).

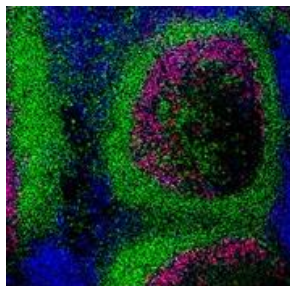

Warranty

The products sold hereunder are warranted only to be free from defects in workmanship and material at the time of delivery to the customer. Miltenyi Biotec GmbH makes no warranty or representation, either expressed or implied, with respect to the fitness of a product for a particular purpose. There are no warranties, expressed or implied, which extend beyond the technical specifications of the products. Miltenyi Biotec GmbH's liability is limited to either replacement of the products or refund of the purchase price. Miltenyi Biotec GmbH is not liable for any property damage, personal injury or economic loss caused by the product.

Miltenyi Biotec GmbH | Friedrich-Ebert-Straße 68 | 51429 Bergisch Gladbach | Germany | Phone +49 2204 8306-0 | Fax +49 2204 85197 | macs@miltenyibiotec.de | www.miltenyibiotec.com
Miltenyi Biotec provides products and services worldwide. Visit www.miltenyibiotec.com/local to find your nearest Miltenyi Biotec contact.

Unless otherwise specifically indicated, Miltenyi Biotec products and services are for research use only and not for therapeutic or diagnostic use. autoMACS, MACS, MACSQuant, Vio, VioBlue, VioBright, and VioGreen are either trademarks or registered trademarks of Miltenyi Biotec GmbH. Copyright © 2016 Miltenyi Biotec GmbH. All rights reserved.