Technical Data Sheet

Alexa Fluor®647 Mouse Anti-Human CD56

Product Information

Material Number: 557711 N-CAM Alternate Name:

100 tests (5µl/test) Size.

5 ul Vol. per Test: B159 Clone:


Mouse IgG1, κ **Isotype:** QC Testing: Human Reactivity:

V NK75 Workshop:

Aqueous buffered solution containing BSA and ≤0.09% sodium azide. Storage Buffer:

Description

Reacts with the 220/135 kDa heavily glycosylated antigen isoforms present on a subpopulation of peripheral blood large granular lymphocytes which demonstrate natural killer activity, but not on myeloid cells, erythrocytes or B cells. This antigen is a pan-NK-cell marker. CD56 is virtually identical to an isoform of the neutral cell adhesion molecule (NCAM), a structure mediating homotypic and heterotypic cell-cell interactions.

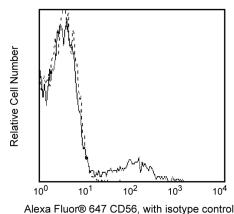

Profile of peripheral blood lymphocytes analyzed by

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number Name Clone Size Alexa Fluor® 647 Mouse IgG1 κ Isotype Control MOPC-21 557714 100 tests

Product Notices

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁶ cells in a 100-μl experimental
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 3. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.

BD Biosciences

bdbiosciences.com

United States Asia Pacific Canada Latin America/Caribbean Europe Japan 877.232.8995 888.259.0187 32.53.720.550 0120.8555.90 65.6861.0633

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



- 5. The Alexa Fluor®, Pacific BlueTM, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific BlueTM dye, and Cascade Blue® dye are covered by pending and issued patents.
- 6. Alexa Fluor is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 8. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 9. Since applications vary, each investigator should titrate the reagent to obtain optimal results.

References

Schlossman S, Boumell L, et al, ed. Leucocyte Typing V. New York: Oxford University Press; 1995.(Biology)
Schlossman SF, Boumsell L, Gilks W, et al, ed. Leukocyte Typing V: White Cell Differentiation Antigens. New York: Oxford University Press; 1995.(Clone-specific)

557711 Rev. 4 Page 2 of 2